



BSI Standards Publication

Sterile urethral catheters for single use

National foreword

This British Standard is the UK implementation of EN ISO 20696:2018. It is identical to ISO 20696:2018. It supersedes BS EN 1616:1997, which is withdrawn.

The UK participation in its preparation was entrusted to Technical Committee CH/84, Catheters and syringes.

A list of organizations represented on this committee can be obtained on request to its secretary.

This publication does not purport to include all the necessary provisions of a contract. Users are responsible for its correct application.

© The British Standards Institution 2020
Published by BSI Standards Limited 2020

ISBN 978 0 539 12735 5

ICS 11.040.25

Compliance with a British Standard cannot confer immunity from legal obligations.

This British Standard was published under the authority of the Standards Policy and Strategy Committee on 31 July 2018.

Amendments/corrigenda issued since publication

Date	Text affected
30 September 2018	Implementation of ISO corrected text September 2018: correction of measurements and units in 5.5.2, 6.4.3, 6.9 and J.2.1; deletion of EN 980 from Bibliography; minor editorial changes
31 January 2020	Implementation of ISO corrected text December 2019: see ISO Foreword for details

English Version

Sterile urethral catheters for single use (ISO 20696:2018)

Sondes urinaires stériles non réutilisables (ISO
20696:2018)

Sterile Harnblasenkatheter zur einmaligen
Verwendung (ISO 20696:2018)

This European Standard was approved by CEN on 5 May 2018.

CEN members are bound to comply with the CEN/CENELEC Internal Regulations which stipulate the conditions for giving this European Standard the status of a national standard without any alteration. Up-to-date lists and bibliographical references concerning such national standards may be obtained on application to the CEN-CENELEC Management Centre or to any CEN member.

This European Standard exists in three official versions (English, French, German). A version in any other language made by translation under the responsibility of a CEN member into its own language and notified to the CEN-CENELEC Management Centre has the same status as the official versions.

CEN members are the national standards bodies of Austria, Belgium, Bulgaria, Croatia, Cyprus, Czech Republic, Denmark, Estonia, Finland, Former Yugoslav Republic of Macedonia, France, Germany, Greece, Hungary, Iceland, Ireland, Italy, Latvia, Lithuania, Luxembourg, Malta, Netherlands, Norway, Poland, Portugal, Romania, Serbia, Slovakia, Slovenia, Spain, Sweden, Switzerland, Turkey and United Kingdom.



EUROPEAN COMMITTEE FOR STANDARDIZATION
COMITÉ EUROPÉEN DE NORMALISATION
EUROPÄISCHES KOMITEE FÜR NORMUNG

CEN-CENELEC Management Centre: Rue de la Science 23, B-1040 Brussels

European foreword

This document (EN ISO 20696:2018) has been prepared by Technical Committee ISO/TC 84 " Devices for administration of medicinal products and catheters " in collaboration with Technical Committee CEN/TC 205 "Non-active medical devices" the secretariat of which is held by DIN.

This European Standard shall be given the status of a national standard, either by publication of an identical text or by endorsement, at the latest by January 2019, and conflicting national standards shall be withdrawn at the latest by January 2019.

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. CEN shall not be held responsible for identifying any or all such patent rights.

This document supersedes EN 1616:1997.

According to the CEN-CENELEC Internal Regulations, the national standards organizations of the following countries are bound to implement this European Standard: Austria, Belgium, Bulgaria, Croatia, Cyprus, Czech Republic, Denmark, Estonia, Finland, Former Yugoslav Republic of Macedonia, France, Germany, Greece, Hungary, Iceland, Ireland, Italy, Latvia, Lithuania, Luxembourg, Malta, Netherlands, Norway, Poland, Portugal, Romania, Serbia, Slovakia, Slovenia, Spain, Sweden, Switzerland, Turkey and the United Kingdom.

Endorsement notice

The text of ISO 20696:2018 has been approved by CEN as EN ISO 20696:2018 without any modification.

Contents

Page

Foreword	v
Introduction	vi
1 Scope	1
2 Normative references	1
3 Terms and definitions	1
4 Intended performance	3
5 General requirements	3
5.1 Risk management	3
5.2 Biocompatibility	3
5.3 Detectability	3
5.4 Surface finish	3
5.5 Size designation	3
5.5.1 General	3
5.5.2 Outer diameter	3
5.5.3 Effective shaft lengths	4
5.6 MRI compatibility	5
5.7 Connector	5
5.8 Sterilization	5
6 Specific requirements	5
6.1 Strength	5
6.2 Connector security	5
6.3 Balloon safety	5
6.4 Catheter inflation lumen integrity and volume maintenance	5
6.4.1 General	5
6.4.2 Compliant balloon	6
6.4.3 Non-compliant balloon	6
6.5 Flow rate	6
6.6 Corrosion resistance	7
6.7 Kink stability	7
6.8 Peak tensile force	7
6.9 Inflated balloon resistance to traction	7
7 Information to be supplied by the manufacturer	8
7.1 General	8
7.2 Marking on the device and/or packaging	8
7.3 Instructions for use	8
Annex A (normative) Test method for determining the strength of the catheter	10
Annex B (normative) Test method for determining the security of fit of the drainage funnel	14
Annex C (normative) Test method for determining balloon safety	16
Annex D (normative) Test method for determining inflation lumen leakage and/or function and/or balloon deflation (catheter with compliant balloon)	18
Annex E (normative) Test method for determination of flow rate through catheter	20
Annex F (normative) Test method for corrosion resistance	22
Annex G (informative) Test method for determining kink stability	23
Annex H (normative) Test method for determining peak tensile force of urethral catheter	25
Annex I (normative) Test method for determining balloon size deflation reliability (catheter with non-compliant balloon)	27
Annex J (normative) Test method for determining inflated balloon resistance to traction	29

Bibliography	33
---------------------------	-----------

Foreword

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

The procedures used to develop this document and those intended for its further maintenance are described in the ISO/IEC Directives, Part 1. In particular the different approval criteria needed for the different types of ISO documents should be noted. This document was drafted in accordance with the editorial rules of the ISO/IEC Directives, Part 2 (see www.iso.org/directives).

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights. Details of any patent rights identified during the development of the document will be in the Introduction and/or on the ISO list of patent declarations received (see www.iso.org/patents).

Any trade name used in this document is information given for the convenience of users and does not constitute an endorsement.

For an explanation on the voluntary nature of standards, the meaning of ISO specific terms and expressions related to conformity assessment, as well as information about ISO's adherence to the World Trade Organization (WTO) principles in the Technical Barriers to Trade (TBT) see the following URL: www.iso.org/iso/foreword.html.

This document was prepared by Technical Committee ISO/TC 84, *Devices for administration of medicinal products and catheters*.

This document is based on EN 1616, *Sterile urethral catheters for single use*.

This corrected version of ISO 20696:2018 incorporates the following corrections:

— in Figure A.1, key items 1 and 2 were inverted.

Introduction

Guidance on transition periods for implementing the requirements of this document is given in ISO/TR 19244.

Sterile urethral catheters for single use

1 Scope

This document specifies requirements and test methods for sterile urethral catheters for single use, with or without a balloon.

This document does not include drainage catheters covered by ISO 20697, e.g. ureteral catheters, nephrostomy catheters, and suprapubic catheters. This document also excludes ureteral stents.

NOTE Ureteral stents are covered in ASTM F1828-97.

2 Normative references

The following documents are referred to in the text in such a way that some or all of their content constitutes requirements of this document. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 14971:—¹⁾, *Medical devices — Application of risk management to medical devices*

ISO 10993-1, *Biological evaluation of medical devices — Part 1: Evaluation and testing within a risk management process*

ISO 80369-1, *Small bore connectors for liquids and gases in healthcare applications — Part 1: General requirements*

3 Terms and definitions

For the purposes of this document, the following terms and definitions apply.

ISO and IEC maintain terminological databases for use in standardization at the following addresses:

- ISO Online browsing platform: available at <https://www.iso.org/obp>
- IEC Electropedia: available at <https://www.electropedia.org/>

3.1

balloon capacity

volume of liquid to be introduced into the catheter in order to fill the inflation channel and inflate the balloon

3.2

coating

substance applied to the surface of the catheter

3.3

compliant balloon

balloon that continues to expand in size as internal pressure increases

3.4

effective length

L_1

length of the catheter that can be inserted into the body

1) To be published (revises ISO 14971:2007). Stage at time of publication ISO/DIS 14971:2013

3.5
effective shaft length

L_3
length of non-perforated portion of the catheter excluding the tip, balloon(s), *funnel(s)* ([3.6](#)), protective sleeves and/or access port(s)

3.6
funnel
proximal portion of the catheter, which may be connected to a drainage system

Note 1 to entry: See [Figures J.1](#) and [J.2](#).

3.7
intermittent
intended to be removed immediately after emptying the bladder

3.8
non-compliant balloon
balloon that expands to one specific size or size range, even as internal pressure increases

3.9
outer diameter
maximum dimension measured across the cylindrical portion of the shaft

3.10
overall length
 L_2
total length from the tip of the catheter to the end of the *funnel* ([3.6](#))

3.11
risk
combination of the probability of occurrence of harm and the severity of that harm

[SOURCE: ISO 14971:—, 3.18]

3.12
risk analysis
systematic use of available information to identify hazards and to estimate the *risk* ([3.11](#))

[SOURCE: ISO 14971:—, 3.19]

3.13
risk assessment
overall process comprising a *risk analysis* ([3.12](#)) and a risk evaluation

[SOURCE: ISO 14971:—, 3.20]

3.14
risk management file
set of records and other documents that are produced by risk management

[SOURCE: ISO 14971:—, 3.25]

3.15
urethral catheter
tubular device intended for being introduced into the urinary bladder through the urethra in order to provide drainage, drug delivery and/or flushing of the bladder

4 Intended performance

The urethral catheter shall demonstrate the ability to accurately and safely access the intended location. The urethral catheter shall demonstrate the ability to drain urine.

5 General requirements

5.1 Risk management

An established risk management process shall be applied to the design of the device.

Compliance shall be checked by inspection of the risk management file verifying compliance to ISO 14971.

5.2 Biocompatibility

The device shall be free from biological hazard in accordance with appropriate testing under ISO 10993-1.

5.3 Detectability

The catheter or at least its effective length shall be detectable by X-ray or by other means (ultra-sound, MRI, etc.), if required by the risk assessment.

NOTE Such as ASTM F640 or DIN 13273-7.

5.4 Surface finish

When examined by normal or corrected to normal vision, the external surface of the effective length of the catheter shall appear free from:

- extraneous matter;
- process and surface defects that may present an unacceptable risk of patient harm.

If deemed necessary based on risk assessment, inspection shall be conducted under a minimum 2,5× magnification.

5.5 Size designation

5.5.1 General

The nominal size of the catheter shall be designated as specified in [5.5.2](#).

5.5.2 Outer diameter

Unless otherwise specified in another clause of this document for a particular type of catheter, the outer diameter shall be expressed as the nominal dimension in millimetres, rounded upwards to the nearest 0,1 mm. Tolerances on this stated size shall be ± 1 French.

NOTE Additional units can also be given. French size (Fr, Ch, FG) is a nominal dimensional identification of the outer size of urethral catheters; calculated as three times the diameter (in millimetres): $Fr = 3 \times D$ (mm).

For devices which are not round by design, the size shall be designated by the dimension of the largest axis. Where relevant, manufacturers may choose to report additional information regarding the device profile, such as the dimension of the second axis for an oval shape.

The balloon capacity shall be expressed in millilitres.

5.5.3 Effective shaft lengths

The minimum effective shaft lengths (L_3) shall be as given in Table 1 (see also Figure 1).

The nominal effective shaft length (L_3) shall be expressed in millimetres, rounded to the nearest millimetre.


NOTE 1 Additional units can also be given.

NOTE 2 Tolerances to the effective length are not specified.

Table 1 — Effective shaft length

Catheter type	L_3 (minimum) mm
Paediatric male	150
Paediatric female	45
Female	60
Male	275

NOTE Shaft dimensions: Shorter sizes can be produced using appropriate risk based clinical justification.

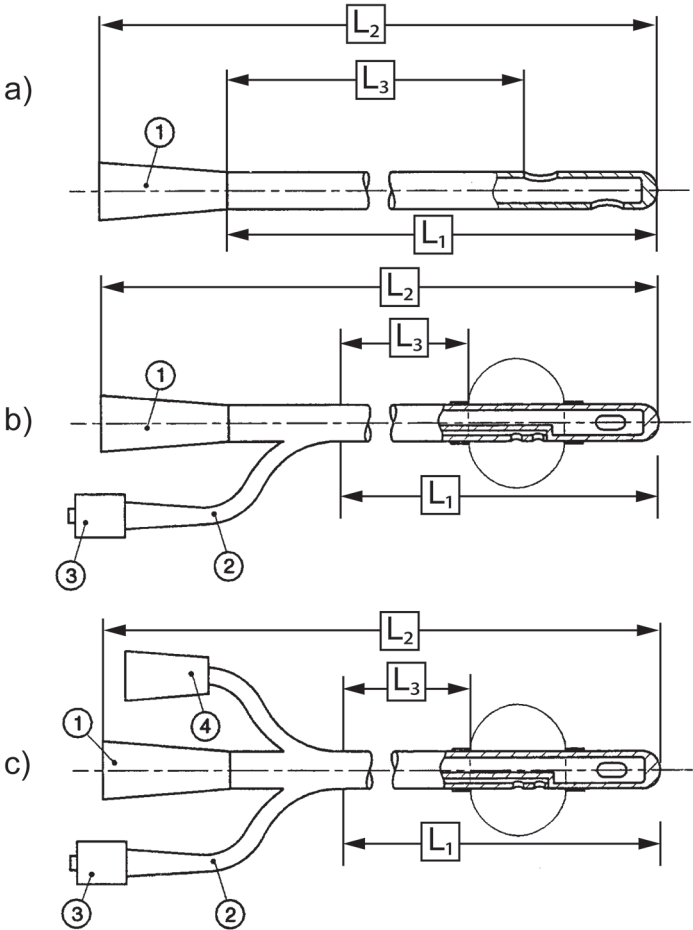

Key

- | | |
|--------------------------|------------------------------|
| 1 drainage funnel | L_1 effective length |
| 2 inflation access port | L_2 overall length |
| 3 valve | L_3 effective shaft length |
| 4 irrigation access port | |

Figure 1 — Typical urethral catheters with and without balloon

5.6 MRI compatibility

If applicable, the hazards of urethral catheters in the magnetic resonance environment shall be evaluated by an appropriate method.

NOTE Such as ASTM F2052, F2213, F2182, and F2119.

5.7 Connector

There is no standardized connector. However, risk of misconnection shall be avoided. This shall be determined by the manufacturer based on risk assessment according to the general requirements of ISO 80369-1.

NOTE 1 A future part of ISO 80369 will address connectors for urethral and urinary applications.

NOTE 2 The funnel is a connecting part, but does not comply with the requirements of ISO 80369-1.

5.8 Sterilization

Urethral catheters and accessories that are sterile shall comply with international, national or regional standards and shall have a sterility assurance level (SAL) of 10^{-6} .

NOTE See applicable parts of ISO 17665, ISO 11135 and ISO 11137 (all parts) for appropriate methods of sterilization.

6 Specific requirements

6.1 Strength

When tested the tip/shaft union and lateral drainage holes shall not show any sign of breaking and neither the tip nor the funnel shall become detached from the shaft.

Compliance shall be checked by the test method in [Annex A](#).

6.2 Connector security

When tested the drainage funnel shall not part from the test connector.

Compliance shall be checked by the test method in [Annex B](#).

6.3 Balloon safety

If present, the balloon shall not leak and shall not occlude the lateral drainage holes.

Compliance shall be checked by the test method in [Annex C](#).

The change in profile at each end of the uninflated balloon should have a smooth transition to the shaft. The balloon should be capable of approximately symmetrical expansion when filled with water at ambient temperature to its specified balloon capacity.

6.4 Catheter inflation lumen integrity and volume maintenance

6.4.1 General

If a balloon is present, choose the appropriate requirement from [6.4.2](#) and [6.4.3](#).