

SAFETY ALERT

CLEANING OF CYLINDERS RETURNED FROM HEALTH CARE FACILITIES DURING A PANDEMIC

As health care providers around the world continue to work to assist patients affected by the respiratory illness Coronavirus Disease 2019 (COVID-19), medical gas containers might become externally contaminated due to their use in health care facilities. It is important that these containers are cleaned in a manner that removes contaminants, does not cause damage to the container, valve, or gas product, and does not lead to a hazardous situation during the refilling and use of the container's contents. The recommendations in this safety alert may be followed until such time that the authority having jurisdiction (for example, World Health Organization [WHO], Public Health Agency of Canada, U.S. Centers for Disease Control and Prevention [CDC]) indicates that the pandemic is no longer a public health threat.

As discussed in CGA P-83, *Guidelines for Cleaning Externally Contaminated Medical Gas Containers*, it is the customer's responsibility to return to the gas supplier medical gas containers that are free from contamination [1]. This responsibility should be communicated and acknowledged by the health care facility.

CGA P-83 addresses cleaning and disinfection procedures for contaminated cylinders and offers cautions against the use of bleach/water or ammonia/water cleaning solutions. Disinfection procedures recommend the use of isopropyl alcohol (IPA) wipes or their equivalent. See Section 5.3 of CGA P-83 for specific information [1].

Likewise, CGA V-12, *Leak Detection Fluids Use with Gas Cylinder Packages* offers cautions against using soap/water or other similar solutions, that could possibly contain ammonia or chlorine as components for leak checking as they can contribute to stress corrosion cracking over a period of time on metallic surfaces [2].

Questions have been recently asked whether compressed gas containers may be washed with soap/water solutions, or disinfected using a dilute bleach/water solution following CDC recommendations if external contamination is suspected when containers are returned from health care facilities or patients' homes *and the customer has not acknowledged their responsibility to return uncontaminated containers or could not adequately clean the container.*

In the United States, the CDC recommends cleaning visibly dirty surfaces with detergent or soap and water and using diluted household bleach solutions, alcohol solutions with at least 70% alcohol, and most common U.S. Environmental Protection Agency (EPA) registered household disinfectants to disinfect.¹

During the COVID-19 crisis or other pandemics, consideration may be given to the use of soap/water solutions to clean, or a dilute bleach/water solution followed by a clean water rinse to disinfect returned containers that may not have been properly cleaned and decontaminated. Such consideration should carefully weigh immediate cleaning concerns against possible longer-term effects on the container and the valve.

When isopropyl alcohol (IPA) or other flammable disinfectants are used as a spray application to disinfect containers, care must be taken to ensure that there is no residual left in the valve outlet prior to the cylinder being refilled.

References

[1] CGA P-83, *Guidelines for Cleaning Externally Contaminated Medical Gas Containers*, Compressed Gas Association, Inc. www.cganet.com

[2] CGA V-12, *Leak Detection Fluids Use with Gas Cylinder Package*, Compressed Gas Association, Inc. www.cganet.com

¹ More information is available from the CDC online: <https://www.cdc.gov/disasters/bleach.html>

PLEASE NOTE:

The information contained in this document was obtained from sources believed to be reliable and is based on technical information and experience currently available from members of the Compressed Gas Association, Inc. and others. However, the Association or its members, jointly or severally, make no guarantee of the results and assume no liability or responsibility in connection with the information or suggestions herein contained. Moreover, it should not be assumed that every acceptable commodity grade, test or safety procedure or method, precaution, equipment or device is contained within, or that abnormal or unusual circumstances may not warrant or suggest further requirements or additional procedure.

This document is subject to periodic review, and users are cautioned to obtain the latest edition. The Association invites comments and suggestions for consideration. In connection with such review, any such comments or suggestions will be fully reviewed by the Association after giving the party, upon request, a reasonable opportunity to be heard. Proposed changes may be submitted via the Internet at our web site, www.cganet.com.

This document should not be confused with federal, state, provincial, or municipal specifications or regulations; insurance requirements; or national safety codes. While the Association recommends reference to or use of this document by government agencies and others, this document is purely voluntary and not binding unless adopted by reference in regulations.

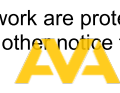
A listing of all publications, audiovisual programs, safety and technical bulletins, and safety posters is available via the Internet at our website at www.cganet.com. For more information contact CGA at Phone: 703-788-2700, ext. 799. E-mail: customerservice@cganet.com.

Work Item 24-55
Cylinder Valve Committee

NOTE—Technical changes from the previous edition are underlined.

SECOND EDITION: 2022
FIRST EDITION: 2020

© 2022 The Compressed Gas Association, Inc. All rights reserved. All materials contained in this work are protected by United States and international copyright laws. You may not alter or remove any trademark, copyright or other notice from this work.

 **Documents**
Document Sharing Hub

CGA GRANTS PERMISSION TO REPRODUCE THIS SAFETY ALERT