



**Health informatics—Personal health device communication**

## **Part 20601: Application profile— Optimized Exchange Protocol**

---

**IEEE Engineering in Medicine and Biology Society**

Sponsored by the  
11073<sup>TM</sup> Standard Committee

---

IEEE  
3 Park Avenue  
New York, NY 10016-5997, USA  
  
19 December 2008

**IEEE Std 11073-20601<sup>TM</sup>-2008**



**Health informatics—Personal health device communication**

# **Part 20601: Application profile— Optimized Exchange Protocol**

Sponsor

**IEEE 11073™ Standard Committee**

of the

**IEEE Engineering in Medicine and Biology Society**

Approved 26 September 2008

**IEEE-SA Standards Board**

**Abstract:** Within the context of the ISO/IEEE 11073 family of standards for device communication, this standard defines a common framework for making an abstract model of personal health data available in transport-independent transfer syntax required to establish logical connections between systems and to provide presentation capabilities and services needed to perform communication tasks. The protocol is optimized to personal health usage requirements and leverages commonly used methods and tools wherever possible.

**Keywords:** medical device communication, personal health devices

The Institute of Electrical and Electronics Engineers, Inc.  
3 Park Avenue, New York, NY 10016-5997, USA

Copyright © 2008 by the Institute of Electrical and Electronics Engineers, Inc.  
All rights reserved. Published 19 December 2008. Printed in the United States of America.

IEEE and IEEE 802 are registered trademarks in the U.S. Patent & Trademark Office, owned by the Institute of Electrical and Electronics Engineers, Incorporated.

**PDF: ISBN 978-0-7381-5826-6 STD95842**  
**Print: ISBN 978-0-7381-5827-3 STDPD95842**

*No part of this publication may be reproduced in any form, in an electronic retrieval system or otherwise, without the prior written permission of the publisher.*

**IEEE Standards** documents are developed within the IEEE Societies and the Standards Coordinating Committees of the IEEE Standards Association (IEEE-SA) Standards Board. The IEEE develops its standards through a consensus development process, approved by the American National Standards Institute, which brings together volunteers representing varied viewpoints and interests to achieve the final product. Volunteers are not necessarily members of the Institute and serve without compensation. While the IEEE administers the process and establishes rules to promote fairness in the consensus development process, the IEEE does not independently evaluate, test, or verify the accuracy of any of the information or the soundness of any judgments contained in its standards.

Use of an IEEE Standard is wholly voluntary. The IEEE disclaims liability for any personal injury, property or other damage, of any nature whatsoever, whether special, indirect, consequential, or compensatory, directly or indirectly resulting from the publication, use of, or reliance upon this, or any other IEEE Standard document.

The IEEE does not warrant or represent the accuracy or content of the material contained herein, and expressly disclaims any express or implied warranty, including any implied warranty of merchantability or fitness for a specific purpose, or that the use of the material contained herein is free from patent infringement. IEEE Standards documents are supplied “AS IS.”

The existence of an IEEE Standard does not imply that there are no other ways to produce, test, measure, purchase, market, or provide other goods and services related to the scope of the IEEE Standard. Furthermore, the viewpoint expressed at the time a standard is approved and issued is subject to change brought about through developments in the state of the art and comments received from users of the standard. Every IEEE Standard is subjected to review at least every five years for revision or reaffirmation. When a document is more than five years old and has not been reaffirmed, it is reasonable to conclude that its contents, although still of some value, do not wholly reflect the present state of the art. Users are cautioned to check to determine that they have the latest edition of any IEEE Standard.

In publishing and making this document available, the IEEE is not suggesting or rendering professional or other services for, or on behalf of, any person or entity. Nor is the IEEE undertaking to perform any duty owed by any other person or entity to another. Any person utilizing this, and any other IEEE Standards document, should rely upon his or her independent judgment in the exercise of reasonable care in any given circumstances or, as appropriate, seek the advice of a competent professional in determining the appropriateness of a given IEEE standard.

**Interpretations:** Occasionally questions may arise regarding the meaning of portions of standards as they relate to specific applications. When the need for interpretations is brought to the attention of IEEE, the Institute will initiate action to prepare appropriate responses. Since IEEE Standards represent a consensus of concerned interests, it is important to ensure that any interpretation has also received the concurrence of a balance of interests. For this reason, IEEE and the members of its societies and Standards Coordinating Committees are not able to provide an instant response to interpretation requests except in those cases where the matter has previously received formal consideration. A statement, written or oral, that is not processed in accordance with the IEEE-SA Standards Board Operations Manual shall not be considered the official position of IEEE or any of its committees and shall not be considered to be, nor be relied upon as, a formal interpretation of the IEEE. At lectures, symposia, seminars, or educational courses, an individual presenting information on IEEE standards shall make it clear that his or her views should be considered the personal views of that individual rather than the formal position, explanation, or interpretation of the IEEE.

Comments for revision of IEEE Standards are welcome from any interested party, regardless of membership affiliation with IEEE. Suggestions for changes in documents should be in the form of a proposed change of text, together with appropriate supporting comments. Comments on standards and requests for interpretations should be submitted to the following address:

Secretary, IEEE-SA Standards Board  
445 Hoes Lane  
Piscataway, NJ 08854  
USA

Authorization to photocopy portions of any individual standard for internal or personal use is granted by The Institute of Electrical and Electronics Engineers, Inc., provided that the appropriate fee is paid to Copyright Clearance Center. To arrange for payment of licensing fee, please contact Copyright Clearance Center, Customer Service, 222 Rosewood Drive, Danvers, MA 01923 USA; +1 978 750 8400. Permission to photocopy portions of any individual standard for educational classroom use can also be obtained through the Copyright Clearance Center.

## Introduction

(This introduction is not part of IEEE Std 11073-20601-2008, Health informatics—Personal health device communication—Part 20601: Application profile—Optimized Exchange Protocol.)

ISO and IEEE 11073 standards enable communication between medical devices and external computer systems. This standard and corresponding IEEE 11073-104zz standards address a need for a simplified and optimized communication approach for personal health devices, which may or may not be regulated devices. These standards align with, and draw upon, the existing clinically focused standards to provide easy management of data from either a clinical or personal health device.

This document addresses a need for an openly defined, independent standard for converting the collected information into an interoperable transmission format so the information can be exchanged between agents and managers.

Other closely related standards include the following:

- ISO/IEEE P11073-00103 [B8]<sup>a</sup> provides an overview of the personal health space and defines the underlying use cases and usage models.
- ISO/IEEE 11073-10101 [B12] documents the nomenclature terms that can be used.
- ISO/IEEE 11073-10201:2004 [B13] documents the extensive domain information model (DIM) leveraged by this standard.
- ISO/IEEE 11073-104zz standards define specific device specializations. For example, ISO/IEEE P11073-10404 [B9] defines how interoperable pulse oximeters work.
- ISO/IEEE 11073-20101:2004 [B14] defines the medical device encoding rules (MDER) used in this standard.

## Notice to users

## Laws and regulations

Users of these documents should consult all applicable laws and regulations. Compliance with the provisions of this standard does not imply compliance to any applicable regulatory requirements. Implementers of the standard are responsible for observing or referring to the applicable regulatory requirements. IEEE does not, by the publication of its standards, intend to urge action that is not in compliance with applicable laws, and these documents may not be construed as doing so.

## Copyrights

This document is copyrighted by ISO/IEEE. It is made available for a wide variety of both public and private uses. These include both use, by reference, in laws and regulations, and use in private self-regulation, standardization, and the promotion of engineering practices and methods. By making this document available for use and adoption by public authorities and private users, ISO/IEEE does not waive any rights in copyright to this document.

---

<sup>a</sup> The numbers in brackets correspond to the numbers of the bibliography in Annex K.

## Updating of IEEE documents

Users of IEEE standards should be aware that these documents may be superseded at any time by the issuance of new editions or may be amended from time to time through the issuance of amendments, corrigenda, or errata. An official IEEE document at any point in time consists of the current edition of the document together with any amendments, corrigenda, or errata then in effect. In order to determine whether a given document is the current edition and whether it has been amended through the issuance of amendments, corrigenda, or errata, visit the IEEE Standards Association website at <http://ieeexplore.ieee.org/xpl/standards.jsp>, or contact the IEEE at the address listed previously.

For more information about the IEEE Standards Association or the IEEE standards development process, visit the IEEE-SA website at <http://standards.ieee.org>.

## Errata

Errata, if any, for this and all other standards can be accessed at the following URL: <http://standards.ieee.org/reading/ieee/updates/errata/index.html>. Users are encouraged to check this URL for errata periodically.

## Interpretations

Current interpretations can be accessed at the following URL: <http://standards.ieee.org/reading/ieee/interp/index.html>.

## Patents

Attention is called to the possibility that implementation of this standard may require use of subject matter covered by patent rights. By publication of this standard no position is taken with respect to the existence or validity of any patent rights in connection therewith. The IEEE is not responsible for identifying Essential Patent Claims for which a license may be required, for conducting inquiries into the legal validity or scope of Patents Claims or determining whether any licensing terms or conditions provided in connection with submission of a Letter of Assurance, if any, or in any licensing agreements are reasonable or nondiscriminatory. Users of this (trial-use, standard, or recommended practice) are expressly advised that determination of the validity of any patent rights, and the risk of infringement of such rights, is entirely their own responsibility. Further information may be obtained from the IEEE Standards Association.

## Participants

At the time this standard was completed, the Personal Health Devices Working Group had the following membership:

**Douglas P. Bogia, *Chair***

Karsten Aalders  
Charles R. Abbruscato  
Maher Abuzaid  
Manfred Aigner  
Murtaza Ali

Deepak Ayyagari  
Merat Bagha  
Doug Baird  
David Baker  
Terry Bartlett

David Bean  
Rudy Belliardi  
Denis Bettini  
Ola Björsne  
Thomas Blackadar

Marc Blanchet  
Douglas P. Bogia  
Terry Bourk  
Bernard Burg  
Lyle G. Bullock, Jr.  
Chris Burns  
Anthony Butt  
Carole C. Carey  
Randy Carroll  
Casper Chen  
James Cheng  
Silviu Chiricescu  
Malcolm Clarke  
Rick A. Cnossen  
Moshe Cohen  
John T. Collins  
Todd Cooper  
Jim DelloStritto  
Matthew d'Entremont  
Kent Dicks  
Jakob Ehrensvar  
Roger M. Ellingson  
Michihiro Enokida  
Mika Erkkilä  
Javier Escayola Calvo  
Leonardo Estevez  
Laurent Falconieri  
Gear Fisher  
Julie N. Fleischer  
Joeseph W. Forler  
Eric Freudenthal  
Miguel Galarraga  
John Garguilo  
Igor Gejdos  
Nicolae Goga  
Chris Gough  
Channa Gowda  
Niclas Granqvist  
Jeff Guttmacher  
Christian Habermann  
Mikey Hagerty  
Rickey L. Hampton  
Sten Hanke  
Jack Harrington  
Kai Hassing  
Hiroshi Hayashi  
Torstein Heggebø  
Ron Hegli  
Rose Higgins  
Kaoru Hiramatsu  
Allen Hobbs  
Kirsten Howard  
Robert Hoy

Robert D. Hughes  
Nick Hunn  
Yutaka Ikeda  
Phil Isaacson  
Ho-In Jeon  
Chris Johnson  
Krishna Jonnalagadda  
Akiyoshi Kabe  
Steve Kahle  
Tomio Kamioka  
Kyung Hee Kang  
Ulf Karlsson  
Andy Kaschl  
Junzo Kashiwara  
Kohichi Kashiwagi  
Ralph Kent  
Kurt Kermes  
Ikuo Keshi  
John Keys  
Alfred Kloos  
Jeongmee Koh  
Alexander Kraus  
Falko Kuester  
Nandu Kushalnagar  
Daniel Lager  
Pierre Landau  
Sungkee Lee  
Yonghee Lee  
Kathryn A. Lesh  
Qiong Li  
Wei-Jung Lo  
Sandra Martinez  
Miguel Martínez de Espronceda  
Cámara  
Jim McCain  
Richard McPartland  
Jinsei Miyazaki  
Florica Moldoveanu  
Brian Møller  
Darr Moore  
Joe Morrissey  
Yoshihiko Motohashi  
Alex Neefus  
Michael E. Nidd  
Hiroaki Niwamoto  
Thomas Norgall  
Yoshiteru Nozoe  
Mikey Paradis  
Jayant Parthasarathy  
Phillip E. Pash  
Thomas Plasa  
Robert E. Ranslam  
Barry Reinhold

Melvin I. Reynolds  
Jeffrey S. Robbins  
Timothy Robertson  
Michael B. Robkin  
Bill Saltzstein  
Stefan Sauermann  
Naveen Saxena  
Paul S. Schluter  
Lars Schmitt  
Mark Schnell  
Richard A. Schrenker  
Aravind Seshagiri  
Marco Sgroi  
Mazen Shihabi  
Robert Smith  
Motoki Sone  
Emily Sopensky  
Ryan Spring  
Nick Steblay  
Costache Stefania  
Lars Steubesand  
John (Ivo) Stivorice  
Ravi Swami  
Xiaorong Tai  
Kunihiro Takiuchi  
Francis Tam  
Haruyuyuki Tatsumi  
Randy Thomas  
Brad Tipler  
Bob Tripp  
Gary Tschautscher  
Masato Tsuchid  
Ken Tubman  
Yoshihiro Uchida  
Sunil Unadkat  
Alpo Värri  
Mark Walters  
Jerry P. Wang  
Jeff Warner  
Toru Watsuji  
Jeff Webber  
Eric White  
David L. Whitlinger  
Vernon C. Williams  
Paul Williamson  
Jan Wittenber  
Ariton Xhafa  
Ricky Yang  
Done-Sik Yoo  
Thomas Zhao  
Daidi Zhong  
Szymon Zysko

The following members of the balloting committee voted on this standard. Balloters may have voted for approval, disapproval, or abstention.

Toru Aihara  
Douglas Bogia  
Lyle Bullock  
Randy Carroll  
Keith Chow  
Malcolm Clarke  
Rick Cnossen  
F. A. Denbrock  
Sourav Dutta  
Julie Fleischer  
Sergiu Goma  
Randall Groves  
Michael Hagerty

Kai Hassing  
Werner Hoelzl  
Robert Hughes  
Philip Isaacson  
Atsushi Ito  
Steven Kahle  
Piotr Karocki  
Kurt Kermes  
Randall Krohn  
Thomas Kurihara  
Qiong Li  
William Lumpkins

Apurva Mody  
Michael S. Newman  
Jayant Parthasarathy  
Melvin Reynolds  
John Rhoads  
Paul Schluter  
Lars Schmitt  
Mark Schnell  
David Simons  
Lars Steubesand  
Eric White  
Jan Wittenber  
Daidi Zhong

When the IEEE-SA Standards Board approved this standard on 26 September 2008, it had the following membership:

**Robert M. Grow**, *Chair*  
**Tom A. Prevost**, *Vice Chair*  
**Steve M. Mills**, *Past Chair*  
**Judith Gorman**, *Secretary*

Victor Berman  
Richard DeBlasio  
Andrew Drozd  
Mark Epstein  
Alexander Gelman  
William R. Goldbach  
Arnold M. Greenspan  
Kenneth S. Hanus

James Hughes  
Richard H. Hulett  
Young Kyun Kim  
Joseph L. Koepfinger\*  
John Kulick  
David J. Law  
Glenn Parsons

Ronald C. Petersen  
Chuck Powers  
Narayanan Ramachandran  
Jon Walter Rosdahl  
Anne-Marie Sahazizian  
Malcolm V. Thaden  
Howard L. Wolfman  
Don Wright

\*Member Emeritus

Also included are the following nonvoting IEEE-SA Standards Board liaisons:

Satish K. Aggarwal, *NRC Representative*  
Michael H. Kelley, *NIST Representative*

Don Messina  
*IEEE Standards Program Manager, Document Development*

Kathryn Cush  
*IEEE Standards Program Manager, Technical Program Development*

# CONTENTS

|                                                                                                  |    |
|--------------------------------------------------------------------------------------------------|----|
| 1. Overview .....                                                                                | 1  |
| 1.1 Scope .....                                                                                  | 1  |
| 1.2 Purpose .....                                                                                | 1  |
| 1.3 Context .....                                                                                | 2  |
| 2. Normative references.....                                                                     | 5  |
| 3. Definitions, acronyms, and abbreviations .....                                                | 5  |
| 3.1 Definitions .....                                                                            | 5  |
| 3.2 Acronyms and abbreviations .....                                                             | 5  |
| 4. Guiding principles .....                                                                      | 6  |
| 5. Introduction to IEEE 11073 personal health devices.....                                       | 7  |
| 5.1 General .....                                                                                | 7  |
| 5.2 Domain information model (DIM) .....                                                         | 8  |
| 5.3 Service model .....                                                                          | 8  |
| 5.4 Communication model.....                                                                     | 8  |
| 6. Personal health device DIM .....                                                              | 8  |
| 6.1 General .....                                                                                | 8  |
| 6.2 Nomenclature usage .....                                                                     | 9  |
| 6.3 Personal health object class definitions .....                                               | 10 |
| 6.3.1 General.....                                                                               | 10 |
| 6.3.2 MDS class .....                                                                            | 12 |
| 6.3.3 Metric class .....                                                                         | 18 |
| 6.3.4 Numeric class.....                                                                         | 23 |
| 6.3.5 RT-SA class .....                                                                          | 26 |
| 6.3.6 Enumeration class .....                                                                    | 27 |
| 6.3.7 PM-store class.....                                                                        | 30 |
| 6.3.8 PM-segment class .....                                                                     | 34 |
| 6.3.9 Scanner classes.....                                                                       | 37 |
| 6.4 Information model extensibility rules.....                                                   | 45 |
| 7. Personal health device service model.....                                                     | 45 |
| 7.1 General .....                                                                                | 45 |
| 7.2 Association service .....                                                                    | 46 |
| 7.3 Object access services.....                                                                  | 46 |
| 7.4 Specific application of object access EVENT REPORT services for personal health devices..... | 47 |
| 7.4.1 General.....                                                                               | 47 |
| 7.4.2 Confirmed and unconfirmed event reports.....                                               | 47 |
| 7.4.3 Configuration event report .....                                                           | 47 |
| 7.4.4 Agent- and manager-initiated measurement data transmission.....                            | 49 |

|                                                               |    |
|---------------------------------------------------------------|----|
| 7.4.5 Variable, fixed, and grouped format event reports ..... | 50 |
| 7.4.6 Single-person and multiple-person event reports .....   | 50 |
| 7.4.7 Temporarily stored measurements .....                   | 51 |
| 8. Communication model .....                                  | 52 |
| 8.1 General .....                                             | 52 |
| 8.2 System context.....                                       | 52 |
| 8.3 Communications characteristics .....                      | 53 |
| 8.3.1 General.....                                            | 53 |
| 8.3.2 Common communications characteristics.....              | 55 |
| 8.3.3 Reliable communications characteristics .....           | 55 |
| 8.3.4 Best-effort communications characteristics .....        | 56 |
| 8.4 State machines .....                                      | 56 |
| 8.4.1 Agent state machine.....                                | 56 |
| 8.4.2 Manager state machine.....                              | 59 |
| 8.4.3 Timeout variables.....                                  | 60 |
| 8.5 Connected procedure .....                                 | 61 |
| 8.5.1 General.....                                            | 61 |
| 8.5.2 Entry conditions .....                                  | 61 |
| 8.5.3 Normal procedures.....                                  | 61 |
| 8.5.4 Exit conditions .....                                   | 61 |
| 8.5.5 Error conditions .....                                  | 62 |
| 8.6 Unassociated procedure .....                              | 62 |
| 8.6.1 General.....                                            | 62 |
| 8.6.2 Entry conditions .....                                  | 62 |
| 8.6.3 Normal procedures.....                                  | 62 |
| 8.6.4 Exit conditions .....                                   | 62 |
| 8.6.5 Error conditions .....                                  | 62 |
| 8.7 Associating procedure .....                               | 63 |
| 8.7.1 General.....                                            | 63 |
| 8.7.2 Entry conditions .....                                  | 63 |
| 8.7.3 Normal procedures.....                                  | 63 |
| 8.7.4 Exit conditions .....                                   | 67 |
| 8.7.5 Error conditions .....                                  | 67 |
| 8.7.6 Test association.....                                   | 67 |
| 8.8 Configuring procedure.....                                | 69 |
| 8.8.1 General.....                                            | 69 |
| 8.8.2 Entry conditions .....                                  | 69 |
| 8.8.3 Normal procedures.....                                  | 69 |
| 8.8.4 Exit conditions .....                                   | 71 |
| 8.8.5 Error conditions .....                                  | 71 |
| 8.9 Operating procedure .....                                 | 72 |
| 8.9.1 General.....                                            | 72 |
| 8.9.2 Entry conditions .....                                  | 72 |
| 8.9.3 Normal procedures.....                                  | 72 |
| 8.9.4 Exit conditions .....                                   | 83 |
| 8.9.5 Error conditions .....                                  | 83 |
| 8.10 Disassociating procedure .....                           | 85 |
| 8.10.1 General.....                                           | 85 |
| 8.10.2 Entry conditions .....                                 | 85 |
| 8.10.3 Normal procedures.....                                 | 85 |
| 8.10.4 Exit conditions .....                                  | 85 |
| 8.10.5 Error conditions .....                                 | 86 |
| 8.11 Message encoding.....                                    | 86 |
| 8.12 Time coordination.....                                   | 86 |

|                                                                   |     |
|-------------------------------------------------------------------|-----|
| 8.12.1 General .....                                              | 86  |
| 8.12.2 Absolute time .....                                        | 86  |
| 8.12.3 Relative time .....                                        | 88  |
| 8.12.4 High-resolution relative time .....                        | 89  |
| 9. Conformance model .....                                        | 89  |
| 9.1 Applicability .....                                           | 89  |
| 9.2 Conformance specification .....                               | 90  |
| 9.3 Implementation conformance statements (ICSs) .....            | 90  |
| 9.4 General conformance .....                                     | 90  |
| 9.4.1 General ICS .....                                           | 91  |
| 9.4.2 Minimum requirements ICS .....                              | 92  |
| 9.4.3 Service support ICS .....                                   | 93  |
| 9.5 Device additions/extensions ICS .....                         | 94  |
| 9.5.1 General additions/extensions ICS .....                      | 94  |
| 9.5.2 Personal health device DIM object and class (POC) ICS ..... | 95  |
| 9.5.3 POC attribute ICS .....                                     | 95  |
| 9.5.4 POC behavior ICS .....                                      | 96  |
| 9.5.5 POC notification ICS .....                                  | 96  |
| 9.5.6 POC nomenclature ICS .....                                  | 97  |
| Annex A (normative) ASN.1 definitions .....                       | 98  |
| Annex B (informative) Scale and range specification example ..... | 130 |
| Annex C (informative) The PM-store concept .....                  | 132 |
| Annex D (informative) Transport profile types .....               | 137 |
| Annex E (normative) State tables .....                            | 140 |
| Annex F (normative) Medical device encoding rules (MDER) .....    | 151 |
| Annex G (informative) Encoded data type definitions .....         | 163 |
| Annex H (informative) Examples .....                              | 182 |
| Annex I (normative) Nomenclature codes .....                      | 190 |
| Annex J (informative) Derivation and modification history .....   | 194 |
| Annex K (informative) Bibliography .....                          | 197 |

## Health informatics—Personal health device communication

# Part 20601: Application profile— Optimized Exchange Protocol

**IMPORTANT NOTICE:** *This standard is not intended to assure safety, security, health, or environmental protection in all circumstances. Implementers of the standard are responsible for determining appropriate safety, security, environmental, and health practices or regulatory requirements.*

*This IEEE document is made available for use subject to important notices and legal disclaimers. These notices and disclaimers appear in all publications containing this document and may be found under the heading “Important Notice” or “Important Notices and Disclaimers Concerning IEEE Documents.” They can also be obtained on request from IEEE or viewed at <http://standards.ieee.org/IPR/disclaimers.html>.*

## 1. Overview

### 1.1 Scope

Within the context of the ISO/IEEE 11073 family of standards for device communication, this standard defines a common framework for making an abstract model of personal health data available in transport-independent transfer syntax required to establish logical connections between systems and to provide presentation capabilities and services needed to perform communication tasks. The protocol is optimized to personal health usage requirements and leverages commonly used methods and tools wherever possible.

### 1.2 Purpose

This document addresses a need for an openly defined, independent standard for converting the information profile into an interoperable transmission format so the information can be exchanged to and from personal telehealth devices and compute engines (e.g., cell phones, personal computers, personal health appliances, and set top boxes).

### 1.3 Context

Figure 1 shows categories and typical devices supporting the personal health space. Agents (e.g., blood pressure monitors, weighing scales, and pedometers) collect information about a person (or persons) and transfer the information to a manager (e.g., cell phone, health appliance, or personal computer) for collection, display, and possible later transmission. The manager may also forward the data to remote support services for further analysis. The information is available from a range of domains including disease management, health and fitness, or aging independently applications.

The communication path between agent and manager is assumed to be a logical point-to-point connection. Generally, an agent communicates with a single manager at any point in time. A manager may communicate with multiple agents simultaneously using separate point-to-point connections.

The overlay shows the focus area of the IEEE 11073™ Personal Health Devices Working Group. The primary concentration is the interface and data exchange between the agents and manager. However, this interface cannot be created in isolation by ignoring the remainder of the solution space. Remaining cognizant of the entire system helps to ensure that data can reasonably move from the agents all the way to the remote support services when necessary. This path may include converting the data format, exchange protocols, and transport protocols across different interfaces. Much of the standardization effort is outside of the scope of the Personal Health Devices Working Group; however, aligning all standardization efforts allows data to flow seamlessly through the overall set of systems.

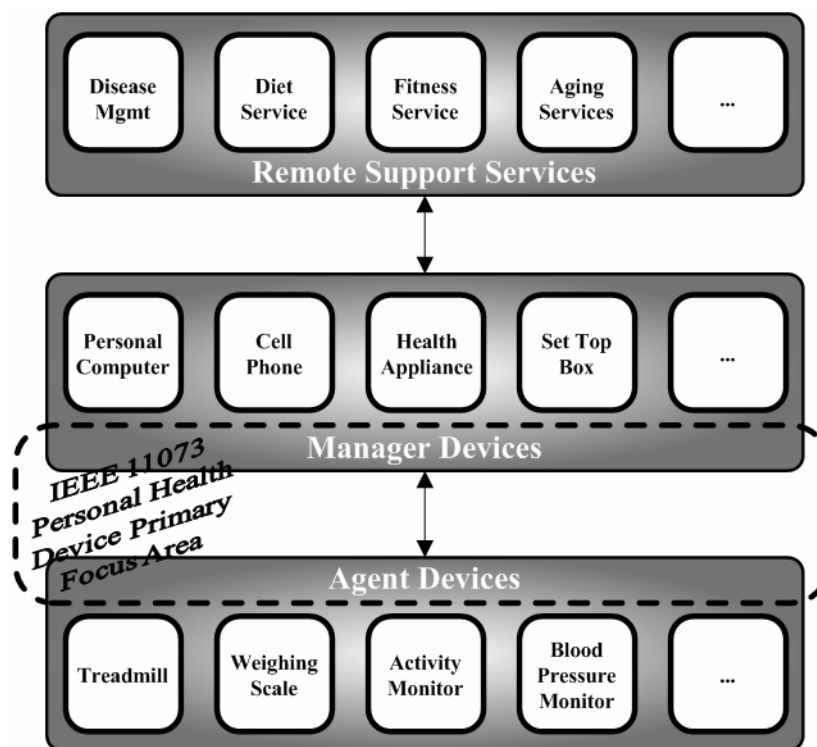


**Figure 1—Overall context of work**

Figure 2 shows a hierarchical view of the architecture of an agent or manager superimposed with a view of the related standards. The application layers are, for the most part, not specific to any particular transport. Where necessary, this standard identifies assumptions that require direct support by a transport or a “shim”

layer above the transport. This approach allows support for various transports. The definition of the transports is outside the scope of this standard and the working group.

Above the transport layer is the Optimized Exchange Protocol (described in this standard). This protocol consists of two aspects: the application layer services and the definition of the data exchange protocol between agents and managers. The application layer services provide the protocol for connection management and reliable transfer of actions and data between agent and manager. The data exchange protocol defines the commands, agent configuration information, data format, and overall protocol. The Optimized Exchange Protocol provides the basis to support any type of agent. For a specific device type, the reader is directed to the device specialization for that agent to understand the capabilities of the device and its implementation according to this standard. The device specialization indicates which aspects of this standard to comprehend and where further information to implement the device is found.

Above the exchange protocol are device specializations that describe specific details relative to the particular agent (e.g., blood pressure monitor, weighing scale, or pedometer). The specializations describe the details of how these agents work and act as a detailed description for creating a specific type of agent. Additionally, they provide reference to a related standard for further details. The standard numbers reserved for device specializations range from IEEE Std 11073-10401 through IEEE Std 11073-10499, inclusive. When the collection of standards is being referenced, the term *IEEE 11073-104zz* is used where *zz* could be any number in the range from 01 to 99, inclusive.

The ISO/IEEE P11073-00103 [B8]<sup>1</sup> technical report describes the overall personal health space with further definition of the underlying use cases and usage models.

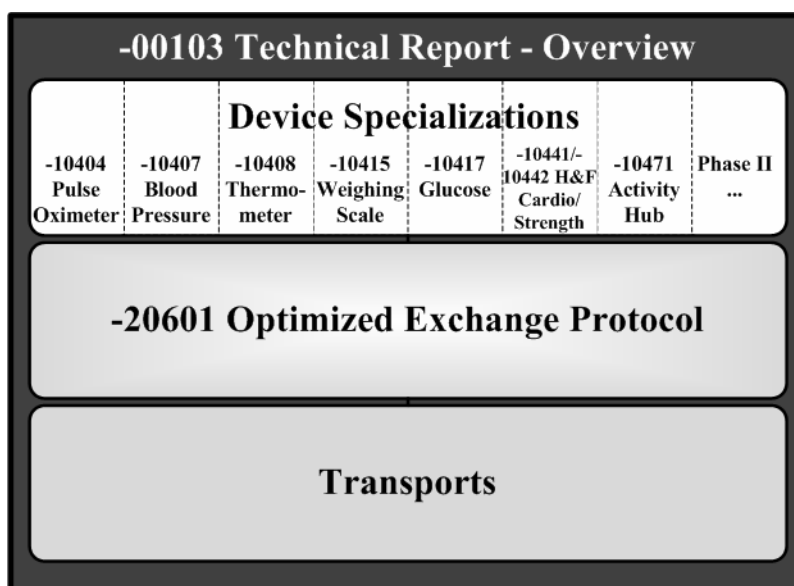


Figure 2—Document map

The personal health device specializations are not being created independently of all other standards. There are a number of existing standards generated for clinical environments upon which these standards draw. Figure 3 shows the relationship to the remainder of the IEEE 11073 documents. There are two types of relationships:

<sup>1</sup> The numbers in brackets correspond to the numbers of the bibliography in Annex K.

- Drawing ideas and/or content from the other documents (dashed lines)
- Leveraging information from the other document and introducing new content into that document to support this standard (solid lines)

This standard imports information from ISO/IEEE 11073-10201:2004 [B13] and ISO/IEEE 11073-20101:2004 [B14] as normative annexes. If there is a discrepancy between these standards, this standard takes priority. Because of the reuse of constructs from these standards, some of the names appear to be more clinically focused [e.g., medical device system (MDS) instead of personal health device system]; however, to maintain consistency, the traditional names have been preserved.

This standard replicates relevant portions of ISO/IEEE 11073-10101 [B12] and incorporates new nomenclature codes.

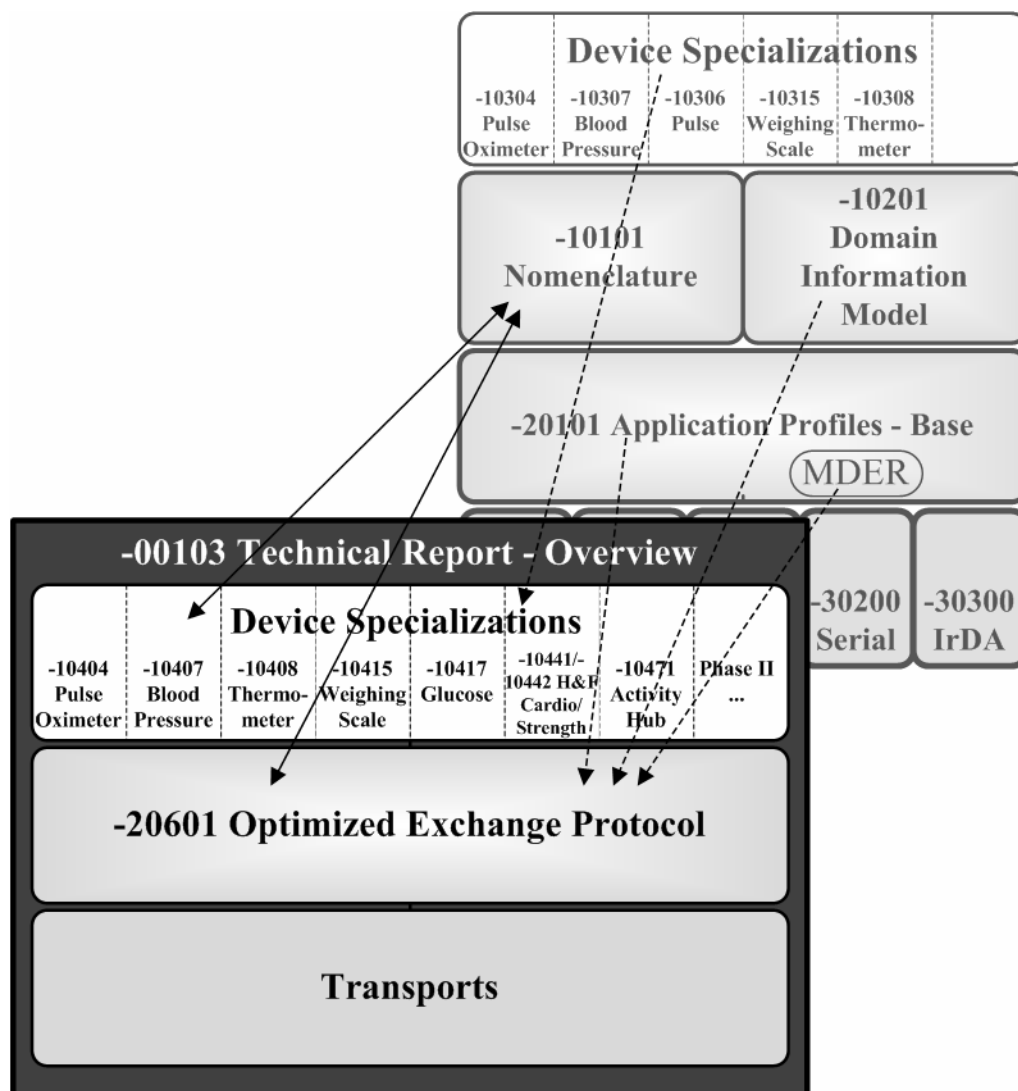


Figure 3—Relationship to other IEEE 11073 documents

## 2. Normative references

The following referenced documents are indispensable for the application of this standard (i.e., they must be understood and used; therefore, each referenced document is cited in the text and its relationship to this standard is explained). For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments or corrigenda) applies.

IEEE Std 802<sup>®</sup>-2001, IEEE Standard for Local and Metropolitan Area Networks: Overview and Architecture.<sup>2</sup>

ITU-T Rec. X.667 (Sept. 2004), Information technology – Open Systems Interconnection – Procedures for the operation of OSI Registration Authorities: Generation and registration of universally unique identifiers (UUIDs) and their use as ASN.1 object identifier components.<sup>3</sup>

## 3. Definitions, acronyms, and abbreviations

### 3.1 Definitions

For the purposes of this standard, the following terms and definitions apply. *The Authoritative Dictionary of IEEE Standards* [B6] should be referenced for terms not defined in this clause.

**3.1.1 agent:** A node that collects and transmits personal health data to an associated manager.

**3.1.2 compute engine:** *See:* **manager**.

**3.1.3 confirmed:** An application-level, completion notification service mechanism. For EVENT REPORT services (i.e., the data plane), confirmation allows the agent to know when the manager has “accepted responsibility” for a piece of data so that the agent can delete that data. For the ACTION, GET, and SET services (i.e., the control plane), confirmation allows the manager to know when the agent has “completed” the requested transaction.

**3.1.4 device:** A physical device implementing either an agent or manager role.

**3.1.5 handle:** An unsigned 16-bit number that is locally unique and identifies one of the object instances within an agent.

**3.1.6 manager:** A node receiving data from one or more agent systems. Examples of managers include a cellular phone, health appliance, set top box, or computer system.

**3.1.7 personal health device:** A device used in personal health applications.

**3.1.8 personal telehealth device:** *See:* **personal health device**.

### 3.2 Acronyms and abbreviations

|       |                                                                 |
|-------|-----------------------------------------------------------------|
| ASCII | American Standard Code for Information Interchange <sup>4</sup> |
| ASN.1 | Abstract Syntax Notation One                                    |
| APDU  | application protocol data unit                                  |

<sup>2</sup> IEEE publications are available from the Institute of Electrical and Electronics Engineers, 445 Hoes Lane, Piscataway, NJ 08854, USA (<http://standards.ieee.org/>).

<sup>3</sup> ITU-T publications are available from the International Telecommunications Union, Place des Nations, CH-1211, Geneva 20, Switzerland/Suisse (<http://www.itu.int/>).

<sup>4</sup> Note that throughout this standard the term ASCII is used to mean the character set as defined in ISO/IEC 646 (1991) [B7].

|                        |                                                                         |
|------------------------|-------------------------------------------------------------------------|
| AVA                    | attribute value assertion                                               |
| BER                    | binary encoding rules                                                   |
| DIM                    | domain information model                                                |
| EUI-64                 | extended unique identifier (64 bits)                                    |
| GMDN                   | Global Medical Device Nomenclature                                      |
| ICS                    | implementation conformance statement                                    |
| ID                     | identifier                                                              |
| LSB                    | least significant bit                                                   |
| MDER                   | medical device encoding rules                                           |
| MDNF                   | medical device numeric format                                           |
| MDS                    | medical device system                                                   |
| MOC                    | medical object class                                                    |
| MSB                    | most significant bit                                                    |
| NaN                    | not a number                                                            |
| NBO                    | network byte order                                                      |
| NRes                   | not at this resolution                                                  |
| NTP                    | Network Time Protocol                                                   |
| OID                    | object identifier                                                       |
| OUI                    | organizationally unique identifier                                      |
| PDU                    | protocol data unit                                                      |
| PER                    | packed encoding rules                                                   |
| POC                    | <i>personal health device domain information model</i> object and class |
| RC <sub>assoc</sub>    | retry count: association procedure                                      |
| RTC                    | real-time clock                                                         |
| RT-SA                  | real-time sample array                                                  |
| SNTP                   | Simple Network Time Protocol                                            |
| TCP                    | Transmission Control Protocol                                           |
| TO <sub>assoc</sub>    | timeout: association procedure                                          |
| TO <sub>ca</sub>       | timeout: confirmed action service                                       |
| TO <sub>cer-mds</sub>  | timeout: confirmed event report service for the MDS object              |
| TO <sub>cer-pms</sub>  | timeout: confirmed event report service for the PM-store object         |
| TO <sub>cer-scan</sub> | timeout: confirmed event report service for the scanner object          |
| TO <sub>clr-pms</sub>  | timeout: confirmed action service to clear the PM-store object          |
| TO <sub>config</sub>   | timeout: configuration procedure                                        |
| TO <sub>cs</sub>       | timeout: confirmed set service                                          |
| TO <sub>get</sub>      | timeout: get service                                                    |
| TO <sub>release</sub>  | timeout: association release procedure                                  |
| TO <sub>sp-mds</sub>   | timeout: special interservice timeout for the MDS object                |
| TO <sub>sp-pms</sub>   | timeout: special segment transfer timeout for the PM-segment object     |
| UTC                    | universal time coordinated                                              |
| UUID                   | universally unique identifier                                           |
| USB                    | universal serial bus                                                    |
| XER                    | Extensible Markup Language (XML) encoding rules                         |

## 4. Guiding principles

This standard and the other personal health device standards fit in the larger context of the ISO/IEEE 11073 family of standards. The full suite of standards enables agents to interconnect and interoperate with managers and with computerized healthcare information systems.

The communication profile defined in this standard takes into account the specific requirements of personal health agents and managers, which are typically used outside a clinical setting, e.g., mobile or in a person's home:

- Personal health agents typically have very limited computing capabilities.

- Personal health agents typically have a fixed configuration, and they are used with a single manager device.
- Personal health agents are frequently battery powered, mobile devices, using a wireless communication link. Therefore, energy efficiency of the protocol is an important aspect.
- Personal health agents are often not permanently active. For example, a weighing scale may provide data only once or twice a day. An efficient connection procedure is needed for minimum overhead for such devices.
- Personal health managers tend to have more processing power, memory, and storage space so the protocol intentionally places more load on the managers.
- Personal health agents and managers convey information that could be useful to clinical professionals. As such, the quality of the data may be considered to have clinical merit even if acquired in a personal health or remote monitoring environment.

The ISO/IEEE 11073 family of standards is based on an object-oriented systems management paradigm. Data (measurement, state, and so on) are modeled in the form of information objects that are accessed and manipulated using an object access service protocol.

To address the unique requirements of personal health devices, a specialized application profile is defined in this standard. This profile leverages concepts from the ISO/IEEE 11073 family of standards and industry best practices to define an optimized communication profile for this domain:

- Where possible, the communication profile is not specific to any particular transport.
- The information model of the communication profile is built on the ISO/IEEE 11073 domain information model (DIM) and includes optimizations where possible.
- An optimized communication protocol is defined to reduce message size, run-time packet construction, and parsing overheads. This is possible due to the lower complexity of the devices in the personal health domain.
- Required definitions for a protocol implementation are included in this standard, rather than referenced. This approach facilitates easier adoption of this standard. In the case of discrepancies between the normative inclusions and a referenced document, this standard takes precedence.

Where possible, versions of this standard are fully backward compatible with at least two major versions.

NOTE—It is expected that any new additions to the DIM or other relevant parts of the ISO/IEEE 11073 family of standards will be adopted and reflected in future revisions of those standards.<sup>5</sup>

## 5. Introduction to IEEE 11073 personal health devices

### 5.1 General

The overall ISO/IEEE 11073 system model is divided into three principal components: the DIM, the service model, and the communications model. These three models work together to represent data, define data access and command methodologies, and communicate the data from an agent to a manager. Because of the tight relationship between the models, the DIM, service model, and communications model are

---

<sup>5</sup> Notes in text, tables, and figures of a standard are given for information only and do not contain requirements needed to implement the standard.

briefly introduced in 5.2, 5.3, and 5.4, respectively, so that when they are described in more detail, in Clause 6, Clause 7, and Clause 8, respectively, the basic concepts are familiar.

## 5.2 Domain information model (DIM)

The DIM, described in detail in Clause 6, characterizes information from an agent as a set of objects. Each object has one or more attributes. Attributes describe measurement data that are communicated to a manager as well as elements that control behavior and report on the status of the agent.

## 5.3 Service model

The service model, described in detail in Clause 7, provides data access primitives that are sent between the agent and manager to exchange data from the DIM. These primitives include commands such as Get, Set, Action, and Event Report.

## 5.4 Communication model

The communication model, described in detail in Clause 8, supports the topology of one or more agents communicating over point-to-point connections to a single manager. For each point-to-point connection, the dynamic system behavior is defined by a connection state machine. The connection state machine defines the states and substates an agent and manager pair passes through, including states related to connection, association, and operation. The communication model also defines in detail the entry, exit, and error conditions for the respective states including various operating procedures for measurement data transmission. The communication model also includes assumptions regarding the underlying communication layers' behavior.

Another function of the communication model is to convert the abstract data modeling (abstract syntax) used in the DIM into a transfer syntax, for example, to binary messages using medical device encoding rules (MDER), that are sent using the communication model.

# 6. Personal health device DIM

## 6.1 General

Personal health devices, within this standard, are defined using an object-oriented model. This DIM defines several classes for modeling an agent. The model describes an agent device as a set of objects that represent the data sources, as the elements that a manager can use to control the behavior of the agent, and as the mechanism the agent uses to report updates to the status of agent representation. Agent device objects have attributes that represent information and status for the object.

Manager devices communicate with agent device objects through the use of well-defined methods, such as GET and SET, and are defined in each subclause describing an object. Information, such as measurements, is sent from agent data objects to the manager device using event reports.

The information model for the domain of personal health devices is an object-oriented model defining data objects of agents, including their attributes and methods. The use of an object-oriented information model supports the following:

- Separation of specification from implementation through the principle of encapsulation
- Support for evolution through the principle of inheritance
- Support for backward compatibility through the principle of polymorphism

The objects derived from classes defined in the information model represent all data that an agent system can communicate to a manager system by means of the application protocol defined in this standard. Such data are modeled in the form of object attributes. Furthermore, the information model defines specific data access services in the form of methods that are used for data exchange between agent and manager systems. These services model the application protocol messages (data access primitives) defined in this standard.

Objects define the structure and the capabilities of the agent system. The manager system accesses these objects to retrieve data and/or to control the agent system. This standard does not define an information model of the manager system.

The information model is a hierarchical model that defines the logical structure and capability of a personal health device. At a top level, the MDS object represents the properties and services of the device itself, independent of its health data capabilities. Properties of the MDS object include attributes for device identification and further technical descriptive and state data. The application-specific data (e.g., health data and measurement data) provided by the personal health device are modeled in the form of further information objects that are logically contained in the MDS object. The set of object attributes, together with this containment relation, describes the configuration and, as such, the capabilities of the personal health device.

Note that while the definitions in this standard make use of object orientation to define the information model, this practice does not imply use of object-oriented technologies (e.g., object-oriented programming languages) for the implementation of this standard in a particular device. The model is used to define data structures and access methods (protocol messages) in a consistent and maintainable way. Conformance to these definitions is at a communication protocol message level only. Specifically, definitions in this standard are optimized to allow very simple agent implementations (e.g., by use of pre-defined transmission templates). Likewise, the implementation of the manager device is free to choose a design that makes use of the information objects versus other design alternatives.

This standard makes use of information classes and objects that are defined in ISO/IEEE 11073-10201:2004 [B13], but adapts these to the domain of personal health device communication in the following ways:

- The definition of attributes that are mandatory, optional, or conditional may be different.
- Additional object services may be defined.
- Additional attributes may be defined.
- Some features of the original model might not be used.

## 6.2 Nomenclature usage

A key aspect of the DIM is that object classes and attributes are referenced using nomenclature codes found in ISO/IEEE 11073-10101 [B12]. By using a consistent nomenclature, interoperability is enhanced as all implementations maintain the same semantic meaning for the numeric codes. Using nomenclature codes also assists with international implementations as the use of strings is reduced.

The ISO/IEEE 11073 nomenclature is defined as a set of context-dependent partitions. The nomenclature code in each context-dependent partition is defined by a 16-bit code that supports up to 65 536 independent

terms per partition. The partitions are referenced by a 16-bit partition code. When the partition of the nomenclature code is defined through context, then it is possible to use only the 16-bit term code. If the context is not defined or a context-independent term code is required, then this situation is specified by a 32-bit code constructed from the 16-bit partition code together with the 16-bit term code. Table 1 shows the partitions that are defined in this standard and/or ISO/IEEE 11073-10101 [B12].

Term codes from 0xF000 – 0xFFFF in each partition in the nomenclature are reserved for private (vendor-specified) nomenclature codes.

For each nomenclature term, ISO/IEEE 11073-10101 [B12] defines a systematic name that explains the term, a unique code value, and a reference identifier (ID). The reference ID has the form MDC\_XXX\_YYY (with MDC referring to “medical device communication”). Throughout this standard, nomenclature terms and nomenclature codes are referenced by the reference ID.

**Table 1—Partitions in the nomenclature**

| Partition number | Nomenclature category                            |
|------------------|--------------------------------------------------|
| 1                | Object-oriented (OO)                             |
| 2                | Supervisory control and data acquisition (SCADA) |
| 3                | Events                                           |
| 4                | Dimensions (units of measurement)                |
| 5                | Virtual attributes                               |
| 6                | Parameter groups                                 |
| 7                | [Body] sites                                     |
| 8                | Infrastructure                                   |
| 9–127            | Reserved                                         |
| 128              | Personal health devices disease management       |
| 129              | Personal health devices health and fitness       |
| 130              | Personal health devices aging independently      |
| 131–254          | Reserved                                         |
| 255              | Return codes                                     |
| 256              | External nomenclature references                 |
| 257–1023         | Reserved                                         |
| 1024             | Private                                          |
| 1025–65 535      | Reserved                                         |

## 6.3 Personal health object class definitions

### 6.3.1 General

The diagram in Figure 4 uses Unified Modeling Language (UML) to represent the information objects of a personal health agent along with class relationships. The top-most object represents the MDS information and its status (see 6.3.2). There are zero or more numeric, real-time sample array (RT-SA), enumeration, scanner, or PM-store objects associated with an MDS object. There are zero or more PM-segments that contain persistent metrics associated with a PM-store. Numeric, RT-SA, and enumeration are derived from a common metric base class that contains common and shared attributes (see 6.3.3). In general, numeric objects represent episodic measurements (see 6.3.4), RT-SA objects represent continuous samples or wave forms (see 6.3.5), enumeration objects represent event annotations (see 6.3.6), and PM-stores (see 6.3.7) along with PM-segments (see 6.3.8) provide a persistent storage mechanism for metrics that are accessed by the manager at a later time. In addition, a scanner object (further defined in 6.3.9) facilitates the reporting of agent-initiated data transfers.

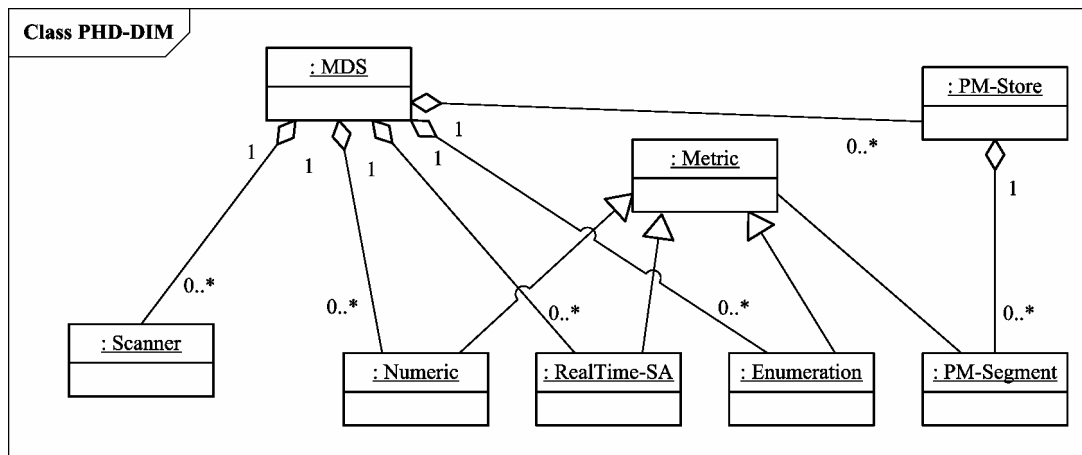


**Figure 4—Personal health device – DIM**

Subclauses 6.3.2 through 6.3.9 describe the classes of the personal health device DIM. Each subclause uses the following format:

- The nomenclature code used to identify the class. This code is used during the configuration event to report the class for each object and allows the manager to learn whether the object being specified is numeric, RT-SA, or any of the other classes.
- The attributes defined by the class.
- The methods available.
- The potential events generated by objects instantiated from the class.
- The available services such as getting or setting attributes.

Each attribute data type is defined using an Abstract Syntax Notation One (ASN.1). The ASN.1 definitions for all data types and exchange formats are found in Annex A.

The attributes for each class are defined in tables that specify the name of the attribute, its nomenclature reference ID, its type, a description of the attribute, and its qualifier. The qualifiers mean O — Attribute is Optional, M — Attribute is Mandatory, and C — Attribute is Conditional and depends on the condition stated in the Remark column. Optional attributes may be implemented on an agent. Mandatory attributes shall be implemented by an agent. Conditional attributes shall be implemented if the condition applies and may be implemented otherwise.

The nomenclature code of the object class (e.g., numeric, RT-SA) is sent to the manager at configuration time to create a mirrored object representation. Each object has a Handle attribute that is used to identify the object for operations (to or from the object) and other attributes to represent and convey information on the physical device and its data sources. Attributes are accessed and modified using methods such as GET and SET. Data are transmitted using EVENTS.

## 6.3.2 MDS class

### 6.3.2.1 General

Each personal health device agent is defined by an object-oriented model as shown in Figure 4. The top-level object of each agent is instantiated from the MDS class. Each agent has one MDS object. The MDS represents the identification and status of the agent through its attributes.

### 6.3.2.2 MDS class identification

The nomenclature code to identify the MDS class is MDC\_MOC\_VMS\_MDS\_SIMP.

### 6.3.2.3 MDS class attributes

Table 2 defines the set of MDS attributes that are supported for personal health agent communication. An MDS object shall support all mandatory attributes, but may have a subset of the conditional and optional attributes present.

**Table 2—MDS attributes**

| Attribute name | Attribute ID       | Attribute type | Remark                                                                                                                                                                                                                                                                                                                                                                                | Qual. |
|----------------|--------------------|----------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Handle         | MDC_ATTR_ID_HANDLE | HANDLE         | The Handle attribute represents a reference ID for this object. The value of the MDS Handle attribute shall be 0.                                                                                                                                                                                                                                                                     | M     |
| System-Type    | MDC_ATTR_SYS_TYPE  | TYPE           | This attribute defines the type of the agent, as defined in nomenclature (e.g., weighing scale). The values shall come from ISO/IEEE 11073-10101 [B12], nom-part-object partition, and subpartition MD-Gen (Medical Device – Generic). Either this attribute or System-Type-Spec-List shall be present. This attribute shall remain unchanged after the configuration is agreed upon. | C     |
| System-Model   | MDC_ATTR_ID_MODEL  | SystemModel    | This attribute defines manufacturer and model number of the agent device. This attribute shall remain unchanged after the configuration is agreed upon.                                                                                                                                                                                                                               | M     |

**Table 2—MDS attributes**

| Attribute name           | Attribute ID               | Attribute type | Remark                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                | Qual. |
|--------------------------|----------------------------|----------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| System-Id                | MDC_ATTR_SYS_ID            | OCTET STRING   | This attribute is an IEEE EUI-64, which consists of a 24-bit unique organizationally unique identifier (OUI) followed by a 40-bit manufacturer-defined ID. The OUI shall be a value assigned by the IEEE Registration Authority ( <a href="http://standards.ieee.org/regauth/index.html">http://standards.ieee.org/regauth/index.html</a> ) and shall be used in accordance with IEEE Std 802-2001. <sup>6</sup> This attribute shall remain unchanged after the configuration is agreed upon.                                                                                        | M     |
| Dev-Configuration-Id     | MDC_ATTR_DEV_CONFIG_ID     | ConfigId       | This attribute defines the identification of the agent device configuration. This Dev-Configuration-Id is static during the lifetime of an association; it is normally exchanged during the association procedure. The manager can GET this attribute during operation. If this attribute is queried prior to when the agent and manager agree upon a configuration, the agent shall return the configuration ID that is being offered at that time. For more information on this attribute, see 8.7.6. This attribute shall remain unchanged after the configuration is agreed upon. | M     |
| Attribute-Value-Map      | MDC_ATTR_ATTRIBUTE_VAL_MAP | AttrValMap     | This attribute defines the attributes that are reported in the fixed format data update messages (see 7.4.5 for more information). Usage of this attribute is mandatory if the agent device uses fixed format value messages to report dynamic data for the object.                                                                                                                                                                                                                                                                                                                   | C     |
| Production-Specification | MDC_ATTR_ID_PROD_SPECN     | ProductionSpec | This attribute defines component revisions, serial numbers, and so on in a manufacturer-specific format. This attribute shall remain unchanged after the configuration is agreed upon.                                                                                                                                                                                                                                                                                                                                                                                                | O     |
| Mds-Time-Info            | MDC_ATTR_MDS_TIME_INFO     | MdsTimeInfo    | This attribute defines the time handling capabilities and the status of the MDS. Usage of this attribute is required if synchronization or settable time is supported.                                                                                                                                                                                                                                                                                                                                                                                                                | C     |

<sup>6</sup> For information on references, see Clause 2.

**Table 2—MDS attributes**

| Attribute name           | Attribute ID              | Attribute type      | Remark                                                                                                                                                                                                                                                                                                                                                                                   | Qual. |
|--------------------------|---------------------------|---------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Date-and-Time            | MDC_ATTR_TIME_ABS         | AbsoluteTime        | This attribute defines the date and time of an agent with resolution of 1/100 of a second, if available. For more information on this attribute, see 8.12. If the agent reports AbsoluteTime in any other message, it shall report its current value of AbsoluteTime in this attribute.                                                                                                  | C     |
| Relative-Time            | MDC_ATTR_TIME_REL         | RelativeTime        | If the agent reports RelativeTime in any other message, it shall report its current value of RelativeTime in this attribute.                                                                                                                                                                                                                                                             | C     |
| HiRes-Relative-Time      | MDC_ATTR_TIME_REL_HI_RES  | HighResRelativeTime | If the agent reports HighResRelativeTime in any other message, it shall report its current value of HighResRelativeTime in this attribute.                                                                                                                                                                                                                                               | C     |
| Date-and-Time-Adjustment | MDC_ATTR_TIME_ABS_ADJUST  | AbsoluteTimeAdjust  | This attribute reports any date and time adjustments that occur either due to a person's changing the clock or events such as daylight savings time. This is used in event reports only. If queried with Get MDS Object command, this value shall be not present or 0. If the agent ever adjusts the date and time, this attribute is used in an event report to report such adjustment. | C     |
| Power-Status             | MDC_ATTR_POWER_STAT       | PowerStatus         | This attribute reports whether power is being drawn from battery or main power lines and the status of charging.                                                                                                                                                                                                                                                                         | O     |
| Battery-Level            | MDC_ATTR_VAL_BATT_CHARGE  | INT-U16             | This attribute reports the percentage of battery capacity remaining, which is undefined if value > 100.                                                                                                                                                                                                                                                                                  | O     |
| Remaining-Battery-Time   | MDC_ATTR_TIME_BATT_REMAIN | BatMeasure          | This attribute represents the predicted amount of operational time left on the batteries. The BatMeasure's unit shall be set to one of MDC_DIM_MIN, MDC_DIM_HR, or MDC_DIM_DAY for minutes, hours, or days, respectively.                                                                                                                                                                | O     |

**Table 2—MDS attributes**

| Attribute name        | Attribute ID                | Attribute type  | Remark                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                          | Qual. |
|-----------------------|-----------------------------|-----------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Reg-Cert-Data-List    | MDC_ATTR_REG_CERT_DATA_LIST | RegCertDataList | This attribute lists various regulatory and/or certification compliance items to which the agent claims adherence as an informative statement. The Implementation Conformance Statements (see Clause 9) take precedence over this attribute and are the legally binding claims. This attribute shall remain unchanged after the configuration is agreed upon.                                                                                                                                                                                                                                                                                                                                                                                   | O     |
| System-Type-Spec-List | MDC_ATTR_SYS_TYPE_SPEC_LIST | TypeVerList     | This attribute reports the type(s) of the agent, as defined in nomenclature (e.g., weighing scale). The values shall come from ISO/IEEE 11073-10101 [B12], nom-part-infrastruct partition, subpartition DEVspec, and reference ISO/IEEE 11073-104zz specializations. If an agent does not follow any specialization, the list shall be left blank. This list shall also contain the version of the specialization. Either this attribute or System-Type shall be present. If the agent is multifunction, this attribute shall be present. This attribute shall remain unchanged after the configuration is agreed upon.                                                                                                                         | C     |
| Confirm-Timeout       | MDC_ATTR_CONFIRM_TIMEOUT    | RelativeTime    | <p>This informational timeout attribute defines the minimum time that the agent shall wait for a Response message from the manager after issuing a Confirmed Event Report invoke message before timing out and transitioning to the Unassociated state.</p> <p>This is an informational attribute for the benefit of the manager. If this attribute is supplied, it shall match the actual timeout value that the agent uses for the Confirmed Event Report generated from the MDS object.</p> <p>This attribute is informational for the manager in the sense that the manager does not use this attribute in an actual implementation of the protocol (i.e., the manager does not time out on an agent-generated Confirmed Event Report).</p> | O     |

**Table 2—MDS attributes**

| Attribute name | Attribute ID | Attribute type | Remark                                                                                                                                             | Qual. |
|----------------|--------------|----------------|----------------------------------------------------------------------------------------------------------------------------------------------------|-------|
|                |              |                | However, the manager might wish to use this information to prioritize its handling of a “short” timeout agent over that of a “long” timeout agent. |       |

The attribute data types are defined in Annex A.

### 6.3.2.4 MDS object methods

Table 3 defines the methods (actions) available for the MDS object. These methods are invoked using the ACTION service. In Table 3, the Method/Action column defines the name of the method. The Mode column defines whether the method is invoked as an unconfirmed action (i.e., roiv-cmip-action from A.10.2) or a confirmed action (i.e., roiv-cmip-confirmed-action). The Action-type column defines the nomenclature ID to use in the action-type field of an action request and response (see A.10.6). The action-info-args column defines the associated data structure to use in the action message for the action-info-args field of the request. The Resulting action-info-args column defines the structure to use in the action-info-args of the response.

**Table 3—MDS object methods**

| Method/Action    | Mode      | Action-type          | action-info-args | Resulting action-info-args |
|------------------|-----------|----------------------|------------------|----------------------------|
| MDS-Data-Request | Confirmed | MDC_ACT_DATA_REQUEST | DataRequest      | DataResponse               |
| Set-Time         | Confirmed | MDC_ACT_SET_TIME     | SetTimeInvoke    | None                       |

— **MDS-Data-Request:**

This method allows the manager system to enable or disable measurement data transmission from the agent (see 8.9.3.3.3 for a description).

— **Set-Time:**

This method allows the manager system to set a real-time clock (RTC) with the absolute time. The agent indicates whether the Set-Time command is valid by using the mds-time-capab-set-clock bit in the Mds-Time-Info attribute (see Table 2).

### 6.3.2.5 MDS object events

Table 4 defines the potential events sent by the MDS object. A manager shall support all methods defined in Table 4.

**Table 4—MDS object events**

| Event                            | Mode                     | Event-type                    | Event-info parameter  | Event-reply-info |
|----------------------------------|--------------------------|-------------------------------|-----------------------|------------------|
| MDS-Configuration-Event          | Confirmed or unconfirmed | MDC_NOTI_CONFIG               | ConfigReport          | ConfigReportRsp  |
| MDS-Dynamic-Data-Update-Var      | Confirmed or unconfirmed | MDC_NOTI_SCAN_REPORT_VAR      | ScanReportInfoVar     | -                |
| MDS-Dynamic-Data-Update-Fixed    | Confirmed or unconfirmed | MDC_NOTI_SCAN_REPORT_FIXED    | ScanReportInfoFixed   | -                |
| MDS-Dynamic-Data-Update-MP-Var   | Confirmed or unconfirmed | MDC_NOTI_SCAN_REPORT_MP_VAR   | ScanReportInfoMPVar   | -                |
| MDS-Dynamic-Data-Update-MP-Fixed | Confirmed or unconfirmed | MDC_NOTI_SCAN_REPORT_MP_FIXED | ScanReportInfoMPFixed | -                |

— **MDS-Configuration-Event:**

This event is sent by the agent during the configuring state of startup if the manager does not already know the agent's configuration from past associations. The event provides static information about the supported measurement capabilities of the agent.

— **MDS-Dynamic-Data-Update-Var:**

This event provides dynamic data (typically measurements) from the agent for some or all of the objects that the agent supports. Data for reported objects are reported using a generic attribute list variable format (see 7.4.5 for details on event report formats). The event is triggered by an MDS-Data-Request from the manager system, or it is sent as an unsolicited message by the agent. For agents that support manager-initiated measurement data transmission, refer to 8.9.3.3.3 for information on controlling the activation and/or period of the data transmission. For agents that do not support manager-initiated measurement data transmission, refer to 8.9.3.3.2 for information on the limited control a manager can assert.

— **MDS-Dynamic-Data-Update-Fixed:**

This event provides dynamic data (typically measurements) from the agent for some or all of the metric objects or the MDS object that the agent supports. Data are reported in the fixed format defined by the Attribute-Value-Map attribute for reported metric objects or the MDS object (see 7.4.5 for details on event report formats). The event is triggered by an MDS-Data-Request from the manager system (i.e., a manager-initiated measurement data transmission), or it is sent as an unsolicited message by the agent (i.e., an agent-initiated measurement data transmission). For agents that support manager-initiated measurement data transmission, refer to 8.9.3.3.3 for information on controlling the activation and/or period of the data transmission. For agents that do not support manager-initiated measurement data transmission, refer to 8.9.3.3.2 for information on the limited control a manager can assert.

— **MDS-Dynamic-Data-Update-MP-Var:**

This is the same as MDS-Dynamic-Data-Update-Var, but allows inclusion of data from multiple persons.

— **MDS-Dynamic-Data-Update-MP-Fixed:**

This is the same as MDS-Dynamic-Data-Update-Fixed, but allows inclusion of data from multiple persons.

### 6.3.2.6 Other MDS services

#### 6.3.2.6.1 GET service

Any agent supporting two-way communication links shall support the GET service to retrieve the values of all implemented MDS object attributes. The GET service can be invoked as soon as the agent receives the

Association Response and moves to the Associated state, including the Operating and Configuring substates.

The manager may request the MDS object attributes of the agent in which case the manager shall send the “Remote Operation Invoke | Get” command (see roiv-cmip-get in A.10.2) with the reserved handle value of 0. The agent shall respond by reporting its MDS object attributes to the manager using the “Remote Operation Response | Get” response (see rors-cmip-get in A.10.2). In the response to a Get MDS Object command, only attributes implemented by the agent are returned.

#### 6.3.2.6.2 SET service

There are currently no settable attributes.

### 6.3.3 Metric class

#### 6.3.3.1 General

The metric class is the base class for all objects representing measurements, status, and context data. The metric class is not instantiated; therefore, it is never part of the agent configuration. As a base class, it defines all attributes, methods, events, and services that are common for all objects representing measurements.

#### 6.3.3.2 Metric class identification

The nomenclature code to identify the metric class is MDC\_MOC\_VMO\_METRIC. This nomenclature code is not used in an agent or a manager implementation as the metric class is just a base class for other classes.

#### 6.3.3.3 Metric class attributes

Table 5 defines the set of metric attributes that are supported for personal health device communication and that are inherited by all metric-derived classes.

**Table 5—Metric attributes**

| Attribute name | Attribute ID       | Attribute type | Remark                                                                                                                                                                                                                                                                         | Qual. |
|----------------|--------------------|----------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Handle         | MDC_ATTR_ID_HANDLE | HANDLE         | The Handle attribute represents a reference ID for this object. Each object shall have a unique ID assigned by the agent. The handle identifies the object in event reports sent to the manager. This attribute shall remain unchanged after the configuration is agreed upon. | M     |

**Table 5—Metric attributes**

| Attribute name         | Attribute ID                 | Attribute type       | Remark                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      | Qual. |
|------------------------|------------------------------|----------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Type                   | MDC_ATTR_ID_TYPE             | TYPE                 | This attribute defines a specific static type of this object as defined in the nomenclature (e.g., pulse rate for a specific numeric object instance). The Type attribute contains the nomenclature partition and term code IDs for context-free, extensible identification. This attribute shall remain unchanged after the configuration is agreed upon.                                                                                                                                                                                                                                                                                                                                  | M     |
| Supplemental-Types     | MDC_ATTR_SUPPLEMENTAL_TYPES  | SupplementalTypeList | This attribute may be used to convey supplemental information about the object beyond the Type and Metric-Id attributes. Supplemental information covers conditions like the location of the sensor or the rate at which the object reacts to changes. Device specializations define the expected usage of this attribute. For example, IEEE Std 11073-10471 [B5] defines location nomenclature for specifying the location of a sensor in a home and ISO/IEEE P11073-10404 [B9] defines three supplemental types for fast response, slow response, and spot checking of the pulse rate or blood oxygenation. This attribute shall remain unchanged after the configuration is agreed upon. | O     |
| Metric-Spec-Small      | MDC_ATTR_METRIC_SPEC_SMALL   | MetricSpecSmall      | This attribute describes the characteristics of the measurements.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                           | M     |
| Metric-Structure-Small | MDC_ATTR_METRIC_STRUCT_SMALL | MetricStructureSmall | This attribute describes the structure of the measurement. If not present, the manager shall assume MetricStructureSmall := {ms-struct-simple, 0}.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                          | O     |
| Measurement-Status     | MDC_ATTR_MSMT_STAT           | MeasurementStatus    | This attribute indicates the validity of a particular value or set of samples.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                              | O     |

**Table 5—Metric attributes**

| Attribute name      | Attribute ID            | Attribute type | Remark                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  | Qual. |
|---------------------|-------------------------|----------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Metric-Id           | MDC_ATTR_ID_PHYSIO      | OID-Type       | <p>This attribute may be used to hold an identification that is more specific than the generic ID in the Type attribute. If the Metric-Id-Partition attribute is valued, it defines the nomenclature partition for this attribute. Otherwise, the OID-Type is taken from the same nomenclature partition as defined in the partition field of the Type attribute.</p> <p>This attribute is needed only if identification changes during operation and the Type attribute does not contain full identification. For example, if the Type attribute contains a generic temperature code (MDC_TEMP), this attribute could report a specific, but changing, identification such as MDC_TEMP_ORAL or MDC_TEMP_RECT. Only one attribute of Metric-Id and Metric-Id-List shall be present.</p> | O     |
| Metric-Id-List      | MDC_ATTR_ID_PHYSIO_LIST | MetricIdList   | <p>This attribute shall be used if a compound observed value is used and does not incorporate the Metric-Id directly (e.g., Compound-Simple-Nu-Observed-Value or Compound-Basic-Nu-Observed-Value) so that elements in the observed value list can be identified individually. The order of the Metric-Id-List shall correspond to the order of the elements in the compound observed value. Only one attribute of Metric-Id and Metric-Id-List shall be present.</p>                                                                                                                                                                                                                                                                                                                   | C     |
| Metric-Id-Partition | MDC_ATTR_METRIC_ID_PART | NomPartition   | <p>This attribute may be used to define the partition from which the Metric-Id or Metric-Id-List nomenclature terms were taken. If not present, the partition is the same as the nomenclature partition defined in the partition field of the Type attribute.</p>                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                       | O     |
| Unit-Code           | MDC_ATTR_UNIT_CODE      | OID-Type       | <p>This attribute defines the nomenclature code for the units of measure from the nom-part-dim partition (e.g., MDC_DIM_KILO_G).</p>                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                    | O     |

**Table 5—Metric attributes**

| Attribute name          | Attribute ID               | Attribute type | Remark                                                                                                                                                                                                                                                                                                                                                                   | Qual. |
|-------------------------|----------------------------|----------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Attribute-Value-Map     | MDC_ATTR_ATTRIBUTE_VAL_MAP | AttrValMap     | This attribute defines the attributes that are reported in the fixed format data update messages. Usage of this attribute is mandatory if the agent uses fixed format value messages to report dynamic measurement data for the object.                                                                                                                                  | C     |
| Source-Handle-Reference | MDC_ATTR_SOURCE_HANDLE_REF | HANDLE         | This attribute establishes a relation of this object instance to a source object (e.g., pulse references sourcing SpO2). This attribute is used whenever it is required to model an explicit relation between object instances to define dependencies. The usage of this attribute is defined by device specializations.                                                 | O     |
| Label-String            | MDC_ATTR_ID_LABEL_STRING   | OCTET STRING   | This attribute defines the textual representation of Type attribute in printable ASCII. The value of this attribute is at the complete discretion of the agent's manufacturer. It could potentially be useful for a manager as a display string or as an aid in deciding how to behave when it does not understand the MDC_ATTR_ID_TYPE as reported by the agent.        | O     |
| Unit-LabelString        | MDC_ATTR_UNIT_LABEL_STRING | OCTET STRING   | This attribute defines the textual representation of Unit-Code dimension in printable ASCII. The value of this attribute is at the complete discretion of the agent's manufacturer. It could potentially be useful for a manager as a display string or as an aid in deciding how to behave when it does not understand the MDC_ATTR_UNIT_CODE as reported by the agent. | O     |

**Table 5—Metric attributes**

| Attribute name        | Attribute ID                   | Attribute type      | Remark                                                                                                                                                                                                                                                                                                                                                                       | Qual. |
|-----------------------|--------------------------------|---------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Absolute-Time-Stamp   | MDC_ATTR_TIME_STAMP_ABS        | AbsoluteTime        | This attribute defines the date and time of observation with resolution of 1/100 of a second, if available. For more information on this attribute, see 8.12. If an agent stores data (either in a PM-store object or as a “temporarily stored measurement”), it shall associate a time stamp (Absolute-Time-Stamp, Relative-Time-Stamp, or HiRes-Time-Stamp) with the data. | C     |
| Relative-Time-Stamp   | MDC_ATTR_TIME_STAMP_REL        | RelativeTime        | This attribute defines the time of observation (time stamp in a relative time format/number of clock ticks as defined by RelativeTime data type). If an agent stores data, it shall associate a time stamp (Absolute-Time-Stamp, Relative-Time-Stamp, or HiRes-Time-Stamp) with the data.                                                                                    | C     |
| HiRes-Time-Stamp      | MDC_ATTR_TIME_STAMP_REL_HI_RES | HighResRelativeTime | This attribute defines the time of observation (time stamp in a high-resolution relative time format/number of clock ticks as defined by a HighResRelativeTime data type). If an agent stores data, it shall associate a time stamp (Absolute-Time-Stamp, Relative-Time-Stamp, or HiRes-Time-Stamp) with the data.                                                           | C     |
| Measure-Active-Period | MDC_ATTR_TIME_PD_MSMT_ACTIVE   | FLOAT-Type          | This attribute defines the time duration of the observation period in seconds.                                                                                                                                                                                                                                                                                               | O     |

#### 6.3.3.4 Metric object methods

None

#### 6.3.3.5 Metric object events

Objects that derive from the metric class do not report their observations directly; rather, the observations are reported through another object, such as the MDS object, a scanner object, or a PM-store object.

### 6.3.3.6 Other metric services

None

### 6.3.4 Numeric class

#### 6.3.4.1 General

An instance of the numeric class represents a numerical measurement. The values of a numeric object are sent from the agent to the manager using the EVENT REPORT service (see 7.3). This class is derived from the metric base class.

#### 6.3.4.2 Numeric class identification

The nomenclature code to identify the numeric class is MDC\_MOC\_VMO\_METRIC\_NU.

#### 6.3.4.3 Numeric class attributes

Table 6 defines the set of numeric attributes that are supported for personal health device communication.

**Table 6—Numeric attributes**

| Attribute name                    | Attribute ID                  | Attribute type      | Remark                                                                                                                                                                                                                                                                                                                                              | Qual. |
|-----------------------------------|-------------------------------|---------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Simple-Nu-Observed-Value          | MDC_ATTR_NU_VAL_OBS_SIMP      | SimpleNuObsValue    | This attribute defines the numerical observed value of the object, without any further embedded status information as found in Nu-Observed-Value. One of Simple-Nu-Observed-Value, Basic-Nu-Observed-Value, Nu-Observed-Value, Compound-Nu-Observed-Value, Compound-Simple-Nu-Observed-Value, or Compound-Basic-Nu-Observed-Value shall be present. | C     |
| Compound-Simple-Nu-Observed-Value | MDC_ATTR_NU_CMPD_VAL_OBS_SIMP | SimpleNuObsValueCmp | This attribute represents an array of Simple-Nu-Observed-Values. One of Simple-Nu-Observed-Value, Basic-Nu-Observed-Value, Nu-Observed-Value, Compound-Nu-Observed-Value, Compound-Simple-Nu-Observed-Value, or Compound-Basic-Nu-Observed-Value shall be present.                                                                                  | C     |

**Table 6—Numeric attributes**

| Attribute name                   | Attribute ID                   | Attribute type     | Remark                                                                                                                                                                                                                                                                                                                                                                                                                     | Qual. |
|----------------------------------|--------------------------------|--------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Basic-Nu-Observed-Value          | MDC_ATTR_NU_VAL_OBS_BASIC      | BasicNuObsValue    | This attribute defines the numerical observed value of the object, without any further embedded status information, but with a smaller numerical representation compared to Simple-Nu-Observed-Value. One of Simple-Nu-Observed-Value, Basic-Nu-Observed-Value, Nu-Observed-Value, Compound-Nu-Observed-Value, Compound-Simple-Nu-Observed-Value, or Compound-Basic-Nu-Observed-Value shall be present.                    | C     |
| Compound-Basic-Nu-Observed-Value | MDC_ATTR_NU_CMPD_VAL_OBS_BASIC | BasicNuObsValueCmp | This attribute represents an array of Basic-Nu-Observed-Values. One of Simple-Nu-Observed-Value, Basic-Nu-Observed-Value, Nu-Observed-Value, Compound-Nu-Observed-Value, Compound-Simple-Nu-Observed-Value, or Compound-Basic-Nu-Observed-Value shall be present.                                                                                                                                                          | C     |
| Nu-Observed-Value                | MDC_ATTR_NU_VAL_OBS            | NuObsValue         | This attribute defines the numerical observed value of the object and combines it with measurement status and unit information. It is used when status/unit are dynamic and are always provided together with the value. One of Simple-Nu-Observed-Value, Basic-Nu-Observed-Value, Nu-Observed-Value, Compound-Nu-Observed-Value, Compound-Simple-Nu-Observed-Value, or Compound-Basic-Nu-Observed-Value shall be present. | C     |

**Table 6—Numeric attributes**

| Attribute name             | Attribute ID             | Attribute type | Remark                                                                                                                                                                                                                                                                                                                                    | Qual. |
|----------------------------|--------------------------|----------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Compound-Nu-Observed-Value | MDC_ATTR_NU_CMPD_VAL_OBS | NuObsValueCmp  | This attribute combines an array of value, status, and unit. This attribute is available for use only in variable format event reports. One of Simple-Nu-Observed-Value, Basic-Nu-Observed-Value, Nu-Observed-Value, Compound-Nu-Observed-Value, Compound-Simple-Nu-Observed-Value, or Compound-Basic-Nu-Observed-Value shall be present. | C     |
| Accuracy                   | MDC_ATTR_NU_ACCUR_MSMT   | FLOAT-Type     | This attribute defines the maximum deviation of actual value from reported observed value (if it can be specified). This attribute shall remain unchanged after the configuration is agreed upon.                                                                                                                                         | O     |

The attributes Compound-Simple-Nu-Obs-Value, Compound-Basic-Nu-Obs-Value, and Compound-Nu-Observed-Value represent a list concept for observed values. This concept should be used whenever a strong relationship between the individual observed values is given, which might be nomenclature and/or application dependent. The compound observed values share the same static attribution context except for the identification of elements. An example is the blood pressure application, where the nomenclature base term expresses “Blood Pressure” and more specific terms express “Blood Pressure Systolic,” “Blood Pressure Diastolic,” and “Blood Pressure Mean.” The corresponding DIM would contain only a single instance of a numeric object, which would use one of the compound numeric observed value formats to represent the “systolic,” “diastolic,” and “mean” parts of the “Blood Pressure.”

#### 6.3.4.4 Numeric object methods

None

#### 6.3.4.5 Numeric object events

None

#### 6.3.4.6 Other numeric services

None

### 6.3.5 RT-SA class

#### 6.3.5.1 General

An instance of the RT-SA class represents a wave form measurement. The values of the RT-SA object are sent from the agent to the manager using the EVENT REPORT service (see 7.3). This class is derived from the metric base class.

#### 6.3.5.2 RT-SA class identification

The nomenclature code to identify the RT-SA class is MDC\_MOC\_VMO\_METRIC\_SA\_RT.

#### 6.3.5.3 RT-SA class attributes

Table 7 defines the set of RT-SA attributes that are supported for personal health device communication.

**Table 7—RT-SA attributes**

| Attribute name                | Attribute ID                                                                    | Attribute type                                          | Remark                                                                                                                                                                                                                                                                                                                                          | Qual. |
|-------------------------------|---------------------------------------------------------------------------------|---------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Sample-Period                 | MDC_ATTR_TIME_PD_SAMP                                                           | RelativeTime                                            | This attributes defines time interval between successive samples given in 1/8 of a millisecond. Thus, 8000 = 1 s. This attribute shall remain unchanged after the configuration is agreed upon.                                                                                                                                                 | M     |
| Simple-Sa-Observed-Value      | MDC_ATTR_SIMP_SA_OBS_VAL                                                        | OCTET STRING                                            | This byte array contains the samples that are reported by the agent in the format that is described by the Sa-Specification and Scale-and-Range-Specification. The length shall be even with padding bytes at the end. Sa-Specification defines the actual number of utilized bytes.                                                            | M     |
| Scale-and-Range-Specification | MDC_ATTR_SCALE_SPECN_I8<br>MDC_ATTR_SCALE_SPECN_I16<br>MDC_ATTR_SCALE_SPECN_I32 | ScaleRangeSpec8<br>ScaleRangeSpec16<br>ScaleRangeSpec32 | This attribute defines mapping between samples and actual values as well as measurement range. The type depends on sample resolution (sample-size field within sample-type field of Sa-Specification). Exactly one of the three specifications shall be included. This attribute shall remain unchanged after the configuration is agreed upon. | M     |
| Sa-Specification              | MDC_ATTR_SA_SPECN                                                               | SaSpec                                                  | This attribute describes the sample array and sample types. This attribute shall remain unchanged after the configuration is agreed upon.                                                                                                                                                                                                       | M     |

Characteristics of the RT-SA object can be gained by examination of the Sa-Specification attribute. This attribute defines the number of elements in the array and the size of an element and is more thoroughly described in A.3.4.

For agents that support manager-initiated measurement data transmission, refer to 8.9.3.3.3 for information on controlling the activation and/or period of the data transmission. For agents that do not support manager-initiated measurement data transmission, refer to 8.9.3.3.2 for information on the limited control a manager can assert.

— **Scale-and-Range-Specification:**

The Scale-and-Range-Specification attribute defines the coefficients for an algorithm to map the scaled values into their absolute values. The manager shall apply the following algorithm:

$$Y = M \times X + B$$

where

Y = the converted absolute value

M = (upper-absolute-value – lower-absolute-value) / (upper-scaled-value – lower-scaled-value)

B = upper-absolute-value – (M × upper-scaled-value)

X = the scaled value

An example of this algorithm in use can be found in Annex B.

Note that the term *absolute-value* does not refer to the mathematical absolute value in which all values are positive, but rather to the actual, measured value.

#### 6.3.5.4 RT-SA object methods

None

#### 6.3.5.5 RT-SA object events

None

#### 6.3.5.6 Other RT-SA services

None

### 6.3.6 Enumeration class

#### 6.3.6.1 General

An instance of the enumeration class represents status information and/or annotation information. The values of the enumeration object are coded in the form of normative codes (as defined in ISO/IEEE 11073-10101 [B12]) or in the form of free text. The values of the enumeration object are sent from the agent to the manager using the EVENT REPORT service (see 7.3). This class is derived from the metric base class.

### 6.3.6.2 Enumeration class identification

The nomenclature code to identify the enumeration class is MDC\_MOC\_VMO\_METRIC\_ENUM.

### 6.3.6.3 Enumeration class attributes

Table 8 defines the set of enumeration attributes that are supported for personal health device communication.

**Table 8—Enumeration attributes**

| Attribute name                     | Attribute ID                        | Attribute type | Remark                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  | Qual. |
|------------------------------------|-------------------------------------|----------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Enum-Observed-Value-Simple-OID     | MDC_ATTR_ENUM_OBS_VAL_SIMP_OID      | OID-Type       | The value is reported as a nomenclature code. If the Enum-Observed-Value-Partition attribute is valued, it defines the nomenclature partition for this attribute. Otherwise, the OID-Type is taken from the same nomenclature partition as defined in the partition field of the Type attribute. One of Enum-Observed-Value-Simple-OID, Enum-Observed-Value-Simple-Bit-Str, Enum-Observed-Value-Basic-Bit-Str, Enum-Observed-Value-Simple-Str, or Enum-Observed-Value shall be present. | C     |
| Enum-Observed-Value-Simple-Bit-Str | MDC_ATTR_ENUM_OBS_VAL_SIMP_BIT_STR  | BITS-32        | The value is reported as a bit string of 32-bits. One of Enum-Observed-Value-Simple-OID, Enum-Observed-Value-Simple-Bit-Str, Enum-Observed-Value-Basic-Bit-Str, Enum-Observed-Value-Simple-Str, or Enum-Observed-Value shall be present.                                                                                                                                                                                                                                                | C     |
| Enum-Observed-Value-Basic-Bit-Str  | MDC_ATTR_ENUM_OBS_VAL_BASIC_BIT_STR | BITS-16        | The value is reported as a bit string of 16-bits. One of Enum-Observed-Value-Simple-OID, Enum-Observed-Value-Simple-Bit-Str, Enum-Observed-Value-Basic-Bit-Str, Enum-Observed-Value-Simple-Str, or Enum-Observed-Value shall be present.                                                                                                                                                                                                                                                | C     |

**Table 8—Enumeration attributes**

| Attribute name                 | Attribute ID                   | Attribute type      | Remark                                                                                                                                                                                                                                                                                                                   | Qual. |
|--------------------------------|--------------------------------|---------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Enum-Observed-Value-Simple-Str | MDC_ATTR_ENUM_OBS_VAL_SIMP_STR | EnumPrintableString | The value is reported as an ASCII printable string. One of Enum-Observed-Value-Simple-OID, Enum-Observed-Value-Simple-Bit-Str, Enum-Observed-Value-Basic-Bit-Str, Enum-Observed-Value-Simple-Str, or Enum-Observed-Value shall be present.                                                                               | C     |
| Enum-Observed-Value            | MDC_ATTR_VAL_ENUM_OBS          | EnumObsValue        | This attribute defines a structured observed value that permits additional flexibility about the data type of the reported value. One of Enum-Observed-Value-Simple-OID, Enum-Observed-Value-Simple-Bit-Str, Enum-Observed-Value-Basic-Bit-Str, Enum-Observed-Value-Simple-Str, or Enum-Observed-Value shall be present. | C     |
| Enum-Observed-Value-Partition  | MDC_ATTR_ENUM_OBS_VAL_PART     | NomPartition        | This attribute may be used to define the partition from which the Enum-Observed-Value-Simple-OID or the Enum-Observed-Value's observation OID nomenclature term was taken. If not present, the partition is the same as the nomenclature partition defined in the partition field of the Type attribute.                 | O     |

#### 6.3.6.4 Enumeration object methods

None

#### 6.3.6.5 Enumeration object events

None

#### 6.3.6.6 Other enumeration services

None

### 6.3.7 PM-store class

#### 6.3.7.1 General

An instance of the PM-store class provides long-term storage capabilities for metric data. Data are stored in a variable number of PM-segment objects (see 6.3.8). The stored data of the PM-store object are requested from the agent by the manager using object access services (see 7.3). Anybody not familiar with the PM-store concept may wish to read Annex C for a conceptual overview prior to reading the following subclauses.

#### 6.3.7.2 PM-store class identification

The nomenclature code to identify the PM-store class is MDC\_MOC\_VMO\_PMSTORE.

#### 6.3.7.3 PM-store class attributes

Table 9 defines the set of PM-store attributes that are supported for personal health device communication:

**Table 9—PM-store attributes**

| Attribute name         | Attribute ID                     | Attribute type | Remark                                                                                                                                                                                                                                                                                                                                                                                                                | Qual. |
|------------------------|----------------------------------|----------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Handle                 | MDC_ATTR_ID_HANDLE               | HANDLE         | The Handle attribute represents a reference ID for this object. Each object shall have a unique ID assigned by the agent. The handle identifies the object in event reports sent to the manager and to address the object instance in messages invoking object methods. This attribute shall remain unchanged after the configuration is agreed upon.                                                                 | M     |
| PM-Store-Capab         | MDC_ATTR_PM_STORE_CAPAB          | PmStoreCapab   | This attribute defines basic capabilities of the PM-store object instance. This attribute shall remain unchanged after the configuration is agreed upon.                                                                                                                                                                                                                                                              | M     |
| Store-Sample-Algorithm | MDC_ATTR_METRIC_STORE_SAMPLE_ALG | StoSampleAlg   | This attribute describes how the sample values stored in the PM-segment have been processed. The StoSampleAlg structure describes the available sampling algorithms. If there is no specific sampling algorithm used (in other words the sample values are raw data), then this attribute shall have a value of st-alg-no-downsampling. This attribute shall remain unchanged after the configuration is agreed upon. | M     |

**Table 9—PM-store attributes**

| Attribute name       | Attribute ID                    | Attribute type   | Remark                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                             | Qual. |
|----------------------|---------------------------------|------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Store-Capacity-Count | MDC_ATTR_METRIC_STORE_CAPAC_CNT | INT-U32          | This attribute is the maximum number of stored PM-segment entries (entries in all contained PM-segments). This attribute shall remain unchanged after the configuration is agreed upon.                                                                                                                                                                                                                                                                                                                                                            | O     |
| Store-Usage-Count    | MDC_ATTR_METRIC_STORE_USAGE_CNT | INT-U32          | This attribute is the actual number of currently stored PM-segment entries (entries in all contained PM-segments).                                                                                                                                                                                                                                                                                                                                                                                                                                 | O     |
| Operational-State    | MDC_ATTR_OP_STAT                | OperationalState | The attribute indicates if new entries are currently being inserted in any of the contained PM-segments. If any PM-segment contained by this PM-store is having data actively added to it, this attribute shall be set to enabled. Otherwise, it shall be set to disabled.                                                                                                                                                                                                                                                                         | M     |
| PM-Store-Label       | MDC_ATTR_PM_STORE_LABEL_STRING  | OCTET STRING     | This attribute is an application-dependent label for the PM-store in printable ASCII to indicate its intended use and may be used for display purposes. This attribute shall remain unchanged after the configuration is agreed upon.                                                                                                                                                                                                                                                                                                              | O     |
| Sample-Period        | MDC_ATTR_TIME_PD_SAMP           | RelativeTime     | This attribute determines the frequency at which entries are added to the PM-segments. This attribute shall be present either in the PM-store (in which case it applies to all periodically storing PM-segments in the PM-store) or alternatively in each PM-segment, if values are sampled periodically, so the time difference for two entries in the Fixed-Segment-Data is constant (i.e., the pmsc-peri-seg-entries bit in the Pm-Store-Capab attribute is set). This attribute shall remain unchanged after the configuration is agreed upon. | C     |
| Number-Of-Segments   | MDC_ATTR_NUM_SEG                | INT-U16          | This attribute is the number of currently instantiated PM-segments contained in the PM-store. Note that the PM-segment attribute Instance-Number is NOT related to this number (i.e., does not need to be in the range from 0 to Number-Of-Segments), but shall be retrieved with the Get-Segment-Info method.                                                                                                                                                                                                                                     | M     |

**Table 9—PM-store attributes**

| Attribute name | Attribute ID           | Attribute type | Remark                                                                                                                                                                                                                                                                                                                                                                                                        | Qual. |
|----------------|------------------------|----------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Clear-Timeout  | MDC_ATTR_CLEAR_TIMEOUT | RelativeTime   | <p>This timeout attribute defines the minimum time that the manager shall wait for the completion of a PM-store clear command.</p> <p>If, after the manager sends a Confirmed Action(Clear Segments) invoke command, the timeout expires before the manager receives the corresponding Confirmed Action response message, the manager shall transition to the Unassociated state as described in 8.9.5.6.</p> | M     |

The attributes Handle and PM-Store-Capab are part of the agent configuration; therefore, the manager knows the corresponding attribute values after the Configuring procedure.

#### 6.3.7.4 PM-store object methods

Table 10 defines the methods (actions) of a PM-store object. These methods can be invoked using the ACTION service.

**Table 10—PM-store object methods**

| Method/Action          | Mode      | Action-type           | action-info-args     | Resulting action-info-args |
|------------------------|-----------|-----------------------|----------------------|----------------------------|
| Clear-Segments         | Confirmed | MDC_ACT_SEG_CLR       | SegmSelection        | (empty)                    |
| Get-Segment-Info       | Confirmed | MDC_ACT_SEG_GET_INFO  | SegmSelection        | SegmentInfoList            |
| Trig-Segment-Data-Xfer | Confirmed | MDC_ACT_SEG_TRIG_XFER | TrigSegmDataXfer-Req | TrigSegmDataXferRsp        |

If an agent supports the PM-store class, the support of the Get-Segment-Info and Trig-Segment-Data-Xfer methods is mandatory. Support for the Clear-Segments method is optional and is indicated in the PM-Store-Capab attribute.

##### — Clear-Segments:

This method allows the manager to delete the data currently stored in one or more selected PM-segments. All entries in the selected PM-segments are deleted. If the agent supports a variable number of PM-segments, the agent may delete empty PM-segments. Additionally, the agent may clear PM-segments without direction from the manager (e.g., the user of the agent could choose to delete data stored on the agent); however, if doing so while in an Associated state, the Instance-Number shall remain empty for the duration of the association. The Instance-Number of all other PM-segments shall be unaffected by clearing a segment. If this method is invoked on a PM-segment that has the Operational-State attribute set to enabled, the agent shall reply with a not-allowed-by-object error (roer) with a return code of MDC\_RET\_CODE\_OBJ\_BUSY.

Note that the behavior of the Clear-Segments method is application specific. The method may remove all entries from the specified PM-segment, leaving it empty, or it may remove the defined PM-segment

completely. This behavior is defined in the PM-Store-Capab attribute. For specific applications, recommendations are defined in corresponding device specializations, making use of the PM-store.

— **Get-Segment-Info:**

This method allows the manager to retrieve PM-segment attributes of one or more PM-segments, with the exception of the Fixed-Segment-Data attribute which contains the actual stored data and is retrieved by using the Trig-Segment-Data-Xfer method. In particular, the Get-Segment-Info method allows the manager to retrieve the Instance-Number attributes of the PM-segment object instances and their data contents.

— **Trig-Segment-Data-Xfer:**

This method allows the manager to start the transfer of the Fixed-Segment-Data attribute of a specified PM-segment. The agent indicates in the response if it accepts or denies this request. If the agent accepts the request, the agent sends Segment-Data-Event messages as described in 6.3.7.5. If this method is invoked on a PM-segment that has the Operational-State attribute set to enabled, the agent shall reply with a not-allowed-by-object error (roer) with a return code of MDC\_RET\_CODE\_OBJ\_BUSY.

### 6.3.7.5 PM-store object events

Table 11 defines the potential events sent by a PM-store object:

**Table 11—PM-store object events**

| Event              | Mode      | Event-type            | Event-info parameter | Event-reply-info  |
|--------------------|-----------|-----------------------|----------------------|-------------------|
| Segment-Data-Event | Confirmed | MDC_NOTI_SEGMENT_DATA | SegmentDataEvent     | SegmentDataResult |

— **Segment-Data-Event:**

This event sends data stored in the Fixed-Segment-Data of a PM-segment from the agent to the manager. The event is triggered by the manager by the Trig-Segment-Data-Xfer method. Once the data transfer is triggered, the agent sends Segment-Data-Event messages until the complete Fixed-Segment-Data is transferred or the transfer is aborted by the manager or agent. See Transfer PM-segment content in 8.9.3.4.2 for a full description.

It is encouraged to place as many segment entries contained in a Segment-Data-Event as possible to reduce the number of messages required for the transfer of the segment.

Support for the event by the agent is mandatory if the agent supports PM-store objects.

### 6.3.7.6 Other PM-store services

#### 6.3.7.6.1 GET service

Support for the GET service shall be provided by any agent that supports one or more PM-store objects only while in the Operating state. The manager uses the GET service to retrieve the values of all PM-store object attributes.

The manager may request the PM-store object attributes of the agent in which case the manager shall send the “Remote Operation Invoke | Get” command (see roiv-cmip-get in A.10.2) with the handle value of the PM-store object, as defined in the agent’s configuration. The agent shall respond by reporting its PM-store

object attributes to the manager using the “Remote Operation Response | Get” response (see rors-cmip-get in A.10.2).

### 6.3.8 PM-segment class

#### 6.3.8.1 General

An instance of the PM-segment class represents a persistently stored episode of measurement data. A PM-segment object is not part of the static agent configuration because the number of instantiated PM-segment instances may dynamically change. The manager accesses PM-segment objects indirectly by methods and events of the PM-store object.

#### 6.3.8.2 PM-segment class identification

The nomenclature code to identify the PM-segment class is MDC\_MOC\_PM\_SEGMENT.

#### 6.3.8.3 PM-segment class attributes

Table 12 defines the set of PM-segment attributes that are supported for personal health device communication.

**Table 12—PM-segment attributes**

| Attribute name       | Attribute ID        | Attribute type    | Remark                                                                                                                                                                                                                                                                                                                                                   | Qual. |
|----------------------|---------------------|-------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Instance-Number      | MDC_ATTR_ID_INSTNO  | InstNumber        | The Instance-Number is the ID of a specific PM-segment object instance. It is used by the manager to address a PM-segment.                                                                                                                                                                                                                               | M     |
| PM-Segment-Entry-Map | MDC_ATTR_PM_SEG_MAP | PmSegmentEntryMap | This attribute defines the format and contents of one stored entry. An entry has an optional header containing information applicable to all elements in the entry. The entry then contains one or more elements, defined by the class, metric ID, handle, and an attribute value map defining the object attributes for each element in the PM-segment. | M     |

**Table 12—PM-segment attributes**

| Attribute name           | Attribute ID                 | Attribute type     | Remark                                                                                                                                                                                                                                                                                                                                                                                                                                             | Qual. |
|--------------------------|------------------------------|--------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| PM-Seg-Person-Id         | MDC_ATTR_PM_SEG_PERSON_ID    | PersonId           | This standard supports devices that have simple support for data from multiple persons. A person ID is used to differentiate different persons. If the PM-store is able to store data for multiple persons, it shall set the pmsc-multi-person bit in the PM-Store-Capab attribute. If this bit is set, all PM-segment instances contained in the PM-store shall support the PM-Seg-Person-Id attribute. Otherwise, this attribute is not defined. | C     |
| Operational-State        | MDC_ATTR_OP_STAT             | OperationalState   | This attribute indicates if new entries are currently being inserted into this PM-segment. If this PM-segment is having data actively added to it, this attribute shall be set to enabled. Otherwise, it shall be set to disabled.                                                                                                                                                                                                                 | M     |
| Sample-Period            | MDC_ATTR_TIME_PD_SAMP        | RelativeTime       | This attribute defines the frequency at which entries are added to the PM-segment. This attribute shall be present either in the PM-store (in which case it applies to all periodically storing PM-segments in the PM-store) or alternatively in each PM-segment. If values are sampled, then the pmsc-peri-seg-entries bit in the PM-Store-Capab attribute shall be set.                                                                          | C     |
| Segment-Label            | MDC_ATTR_PM_SEG_LABEL_STRING | OCTET STRING       | This attribute is an application-dependent label in printable ASCII for the segment to indicate its intended use and may be used for display purposes.                                                                                                                                                                                                                                                                                             | O     |
| Segment-Start-Abs-Time   | MDC_ATTR_TIME_START_SEG      | AbsoluteTime       | This attribute defines the start time of segment.                                                                                                                                                                                                                                                                                                                                                                                                  | O     |
| Segment-End-Abs-Time     | MDC_ATTR_TIME_END_SEG        | AbsoluteTime       | This attribute defines the end time of segment.                                                                                                                                                                                                                                                                                                                                                                                                    | O     |
| Date-and-Time-Adjustment | MDC_ATTR_TIME_ABS_ADJUST     | AbsoluteTimeAdjust | This attribute reports any date and time adjustments that occur either due to a person's changing the clock or events such as daylight savings time. If the agent ever adjusts the date and time, this attribute is reports such adjustment.                                                                                                                                                                                                       | C     |

**Table 12—PM-segment attributes**

| Attribute name      | Attribute ID             | Attribute type                                                                                                           | Remark                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  | Qual. |
|---------------------|--------------------------|--------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Segment-Usage-Count | MDC_ATTR_SEG_USAGE_CNT   | INT-U32                                                                                                                  | This attribute gives the actual (current) number of stored entries.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     | O     |
| Segment-Statistics  | MDC_ATTR_SEG_STATS       | SegmentStatistics                                                                                                        | This attribute defines the array for reporting minimum, mean, maximum statistics for each element to be tagged.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                         | O     |
| Fixed-Segment-Data  | MDC_ATTR_SEG_FIXED_DATA  | N/A<br>The data is stored internal to the device and so this data type never occurs in any protocol definition directly. | This attribute defines the segment data transferred as an array of entries in a format as specified in the PM-Segment-Entry-Map attribute. This is defined here as an opaque data structure without a defined data type. Note that this attribute is not directly accessible; it is only retrievable by the manager using the PM-store Trig-Segm-Data-Xfer method.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      | M     |
| Confirm-Timeout     | MDC_ATTR_CONFIRM_TIMEOUT | RelativeTime                                                                                                             | <p>This informational timeout attribute defines the minimum time that the agent shall wait for a Response message from the manager after issuing a Confirmed Event Report invoke message before timing out and transitioning to the Unassociated state.</p> <p>This is an informational attribute for the benefit of the manager. If this attribute is supplied, it shall match the actual timeout value that the agent uses for the Confirmed Event Report generated from the PM-store object.</p> <p>This attribute is informational for the manager in the sense that the manager does not use this attribute in an actual implementation of the protocol (i.e., the manager does not time out on an agent-generated Confirmed Event Report). However, the manager might wish to use this information to prioritize its handling of a “short” timeout agent over that of a “long” timeout agent.</p> | O     |

**Table 12—PM-segment attributes**

| Attribute name   | Attribute ID              | Attribute type | Remark                                                                                                                                                                                                                                                                                                   | Qual. |
|------------------|---------------------------|----------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|
| Transfer-Timeout | MDC_ATTR_TRANSFER_TIMEOUT | RelativeTime   | <p>This timeout attribute defines the minimum time that the manager shall wait for the complete transfer of PM-segment information.</p> <p>If the timeout expires prior to the reception of the complete PM-segment, the manager shall transition to the Unassociated state as described in 8.9.5.6.</p> | M     |

The Fixed-Segment-Data attribute stores an array of identically formatted entries. In cases where a measurement is not available at the required time, then the value for a numeric measurement represented by the (S)FLOAT-type data type shall use the special NaN (not a number) value to indicate an unavailable value.

The Fixed-Segment-Data attribute may hold very large amounts of data, depending on the agent capabilities and the application. An agent may choose to restrict the maximum size of the Fixed-Segment-Data attribute in a way that is aligned with the maximum transmission unit of the transport system. In order to support this type of behavior, a manager shall be able to support the transfer of Fixed-Segment-Data attributes in multiple application messages.

#### 6.3.8.4 PM-segment object methods

None

#### 6.3.8.5 PM-segment object events

None

#### 6.3.8.6 Other PM-segment services

None

### 6.3.9 Scanner classes

#### 6.3.9.1 General

A scanner object is an observer and “summarizer” of object attribute values. It observes attributes of metric objects (e.g., numeric objects) and generates summaries in the form of notification event reports. See Figure 5 for the class hierarchy of the scanner classes. Each class is described in 6.3.9.2 through 6.3.9.5, respectively.

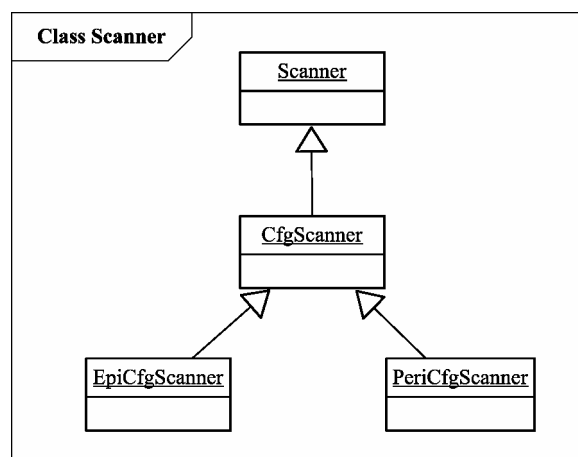


**Figure 5—Personal health device – DIM – scanner model**

The different scanner classes (periodic and episodic) as well as the different instances should be used to distribute different data types across one or many data flows represented by a scanner instance. A pulse oximetry application might choose to have a periodic configurable scanner for RT-SA objects with a Reporting-Interval of 50 ms, a periodic configurable scanner with a Reporting-Interval of 1 s for numeric and enumeration objects, and an episodic configurable scanner for beat-to-beat metric objects (numeric or enumeration objects).

### 6.3.9.2 Scanner class

#### 6.3.9.2.1 General

The scanner class is an abstract class defining attributes, methods, events, and services that are common for its subclasses. As such, it cannot be instantiated.

The scanner concept provides three different event report notifications: variable format, fixed format, and grouped format. See 7.4.5 for the reporting of observed object attributes. The event report formats are described further in 6.3.9.4.5, 6.3.9.5.5, and A.11.5, respectively.

More specialized scanner classes are derived from the scanner base class.

#### 6.3.9.2.2 Scanner class identification

The nomenclature code to identify the scanner class is MDC\_MOC\_SCAN.

#### 6.3.9.2.3 Scanner class attributes

Table 13 defines the set of scanner attributes that are supported for personal health device communication.