

ISO

INTERNATIONAL ORGANIZATION FOR STANDARDIZATION

ISO RECOMMENDATION R 760

DETERMINATION OF WATER BY THE KARL FISCHER METHOD

1st EDITION
June 1968

COPYRIGHT RESERVED

The copyright of ISO Recommendations and ISO Standards belongs to ISO Member Bodies. Reproduction of these documents, in any country, may be authorized therefore only by the national standards organization of that country, being a member of ISO.

For each individual country the only valid standard is the national standard of that country.

Printed in Switzerland

Also issued in French and Russian. Copies to be obtained through the national standards organizations.

BRIEF HISTORY

The ISO Recommendation R 760, *Determination of water by the Karl Fischer method*, was drawn up by Technical Committee ISO/TC 47, *Chemistry*, the Secretariat of which is held by the Ente Nazionale Italiano di Unificazione (UNI).

Work on this question led, in 1962, to the adoption of a Draft ISO Recommendation.

In November 1963, this Draft ISO Recommendation (No. 660) was circulated to all the ISO Member Bodies for enquiry. It was approved, subject to a few modifications of an editorial nature, by the following Member Bodies :

Australia	Hungary	Portugal
Austria	India	Romania
Belgium	Israel	Spain
Chile	Italy	U.A.R.
Colombia	Japan	United Kingdom
Czechoslovakia	Korea, Rep. of	U.S.S.R.
France	Netherlands	Yugoslavia
Germany	Poland	

Two Member Bodies opposed the approval of the Draft :

New Zealand
U.S.A.

The Draft ISO Recommendation was then submitted by correspondence to the ISO Council which decided, in June 1968, to accept it as an ISO RECOMMENDATION.

CONTENTS

	Page
1. Scope	3
2. Principle	3
3. Reactions	3
4. Reagents and auxiliary materials	3
5. Apparatus	5
6. Visual method	6
7. Direct electrometric titration	8
8. Electrometric back titration	9
9. Test report	11
ANNEX A	13
ANNEX B	14
ANNEX C	15

DETERMINATION OF WATER BY THE KARL FISCHER METHOD *

1. SCOPE

This ISO Recommendation describes a method suitable for the determination of free water or water of crystallisation in most solid or liquid chemical products, both organic and inorganic.

Precautions are necessary in certain cases and these are specified in the appropriate ISO Recommendations.


There are two methods of titration, depending on whether the end point is detected visually or electrometrically. The visual method can be used only with colourless solutions and when no electrometric apparatus is available; it is always a direct titration. The electrometric method on the other hand may involve either a direct titration or a back titration. The electrometric method, whether by direct titration or back titration, is the more accurate and for this reason it is recommended.

2. PRINCIPLE

Reaction of any water present with a solution of iodine and sulphur dioxide in a pyridine methanol mixture (Karl Fischer reagent), previously standardized by titration with an exactly known mass of water (see clauses 6.1, 7.1 and 8.1).

NOTE. — Methanol may be replaced by 2-methoxyethanol (ethylene glycol monomethyl ether). With this solvent a more constant titration volume is obtained and the reagent can be used with aldehydes and ketones, without using any special technique.**

3. REACTIONS***

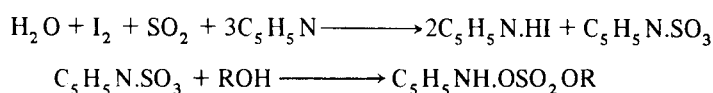

4. REAGENTS AND AUXILIARY MATERIALS

- 4.1 *Methanol*, analytical reagent grade not containing more than 0.05 % of water. If the reagent contains more than this quantity of water, dry it by distillation from magnesium turnings activated with iodine. Collect the distillate in a receiver protected from atmospheric moisture by means of a guard tube filled with anhydrous aluminium sodium silicate or activated silica gel (4.9).
- 4.2 *2-Methoxyethanol (ethylene glycol monomethyl ether)*, not containing more than 0.05 % of water. If the reagent contains more than this quantity of water, dry it by distillation, rejecting the first portion of distillate which will contain the water.

* Fischer, Karl. *Angewandte Chemie* 1935, 48, page 394.

** Peters, E.D. and Jungnickel, I.L. *Anal. Chem.* 1955, 27, 450-53.

*** Smith, D.M. Bryant, W.N.D. and Mitchell Jr., J. *Journal of the American Chemical Society* 1939, 61, page 2407.

- 4.3 *Pyridine*, not containing more than 0.05 % of water. If the reagent contains more than this quantity of water, dry it by distillation, rejecting the first portion of distillate which will contain the water.
- 4.4 *Sample solvent*. Either methanol (4.1), a mixture containing 4 parts of methanol (4.1) and 1 part of pyridine (4.3) (by volume) or (preferably for determinations with compounds containing carbonyl groups) a mixture containing 4 parts of 2-methoxyethanol (4.2) and 1 part of pyridine (4.3) (by volume). In special cases, other solvents may be recommended, for example, acetic acid, pyridine or a mixture containing 1 part of methanol (4.1) and 3 parts of chloroform (by volume).

- 4.5 *Karl Fischer reagent*. Place 670 ml of methanol (4.1) or 2-methoxyethanol (4.2) in a previously dried flask, coloured brown or painted black on the outside, fitted with a ground glass stopper and having a capacity slightly greater than 1 litre.

Add about 85 g of iodine. Stopper the flask and shake it occasionally until the iodine is completely dissolved. Then add approximately 270 ml of pyridine (4.3), stopper the flask again and mix thoroughly. Dissolve 65 g of sulphur dioxide in this solution, cooling to ensure that the temperature of the liquid does not exceed 20 °C.

NOTE. — The reaction being exothermic, it is necessary to cool the flask from the beginning and to maintain it at about 0 °C, for example, by immersing in an ice bath or in crushed solid carbon dioxide.

Replace the ground glass stopper by an attachment for introducing sulphur dioxide. This consists of a cork with a thermometer and an inlet glass tube 6 mm × 8 mm, reaching to within 10 mm of the bottom of the flask, and a small capillary tube for connecting to the atmosphere.

Place the whole assembly with the ice bath on a balance and weigh to the nearest 1 g. Connect the inlet tube to a siphon of sulphur dioxide by means of a flexible connection and a drying tube filled with anhydrous aluminium sodium silicate (4.9) as desiccant and gently open the tap on the siphon.

Adjust the rate of flow of sulphur dioxide so that all the gas is absorbed without the liquid showing any sign of rising in the inlet tube.

Maintain the equilibrium of the balance by gradually increasing the tare and ensure that the temperature of the liquid does not rise above 20 °C. Close the tap on the siphon as soon as the increase in mass reaches 65 g.

Immediately remove the flexible connection and reweigh the flask and its inlet attachment. The mass of dissolved sulphur dioxide should be between 60 and 70 g. A slight excess is not harmful.

Stopper the flask, mix the solution and leave for at least 24 hours before using it. In fact, as the result of imperfectly understood reactions which occur in the fresh reagent, the water equivalent of the reagent decreases rapidly to begin with and then much more slowly.

This water equivalent is between 3.5 and 4.5 mg/ml. It should be determined daily if methanol has been used, but may be done less frequently if 2-methoxyethanol has been used.

It is possible to prepare the Karl Fischer reagent with a lower water content by diluting the solution prepared as described above with the sample solvent (4.4).

Store the reagent out of the light and protected from atmospheric moisture. It should preferably be stored in a reagent bottle of brown or black-painted glass.

- 4.6 *Sodium tartrate*, crystalline, ($\text{Na}_2\text{C}_4\text{H}_4\text{O}_6 \cdot 2\text{H}_2\text{O}$) or water.

- 4.7 *Water/methanol*, 10 mg/ml standard solution.

Using a microburette or a pipette, place 1 ml of water in a dry 100 ml one-mark volumetric flask, containing approximately 50 ml of methanol (4.1). Dilute to the mark with methanol (4.1) and mix (for standardization with this solution see Annex A, clause A.1 or A.2 depending on whether the visual method or the direct electrometric method is used).

4.8 *Water/methanol* solution, approximately 2 g/l.

Using a microburette or a pipette, place 1 ml of water into a perfectly dry 500 ml one-mark volumetric flask containing approximately 100 ml of methanol (4.1). Dilute to the mark with methanol (4.1) and mix (see the correspondence by volume of this solution with the Karl Fischer reagent in clause 8.2.2).

4.9 *Aluminium sodium silicate*, anhydrous, granules of diameter of 1.7 mm for use as a desiccant. These granules may be regenerated by washing with water and drying at 350 °C for at least 48 hours.

Or *activated silica gel*.

4.10 *Silicone base grease*, for lubricating the ground glass joints.

5. APPARATUS

5.1 For direct titration (by visual method or electrometric method).

5.1.1 A suitable form of apparatus for this method is given in Annex B, should no commercial apparatus be available. This standard apparatus comprises the parts described below.

5.1.1.1 *Automatic burette*, 25 ml capacity with a fine point, graduated in 0.05 ml and protected from atmospheric moisture by a guard tube filled with the desiccant (4.9).

5.1.1.2 *Titration vessel*, effective capacity of 100 ml, connected to the tap of the burette (5.1.1.1) by a ground glass joint and having two side tubes, one permitting the introduction of platinum electrodes if the electrometric method is used and the other fitted with a "vaccine cap" to permit the introduction of liquid test samples by means of a syringe without opening the vessel.

5.1.1.3 *Platinum electrodes* (Fig. 1 and 2, Annex B) welded to a glass tube enabling them to be introduced into the bottom of the titration vessel (5.1.1.2) and joined to two copper wires which in turn connect to the device for the electrometric detections of the end point (5.1.1.7) (these are omitted in the case of the visual method).

5.1.1.4 *Electromagnetic stirrer*, operating at 150 to 300 revolutions per minute with a mild steel bar coated with glass or PTFE, fixed on a base of adjustable height.

5.1.1.5 *Reagent bottle for the Karl Fischer reagent* (capacity approximately 3 litres) of brown or black-painted glass in which the filling tube of the automatic burette is immersed through the ground glass stopper.

5.1.1.6 *Rubber bulb* connected to an air dryer filled with the desiccant (4.9) for the admission of dry air under pressure into the reagent bottle (5.1.1.5) in order to fill the burette (5.1.1.1).

5.1.1.7 *Device* for the electrometric detection of the end point shown diagrammatically in Figure 3, Annex B (to be omitted in the case of the visual method).

5.1.2 *Medical syringes*, of suitable capacity and calibrated.

5.1.3 *A small glass tube*, closed at one end and fitted at the other with a rubber stopper, used for weighing and introducing into the titration vessel, for example, the crystalline sodium tartrate (approximately 250 mg) used to standardize the Karl Fischer reagent, or possibly test samples of solid products.

NOTE. — All glassware used should previously be dried at 130 °C in an oven, then placed in a desiccator containing the desiccant (4.9).