

KES

**Analysis of Component
for Diesel Fuel Oil**

KES T – M035

(EST.) 1988 – 07 – 02

Contents

- 1. Scope of Application**
- 2. Gas Chromatograph A Method**
- 3. Gas Chromatograph B Method**



1. Scope of application

This standard shall be prescribed to the method of qualitative and quantitative analysis of an organic component.

2. Gas chromatograph test

2.1 Condition of analysis

- (1) Column : OV-101 0.25 mm ϕ $\times$ 28 m Glass WCOT
- (2) Temperature of column : 70 $\rightarrow$ 230 $^{\circ}\text{C}$ 4 $^{\circ}\text{C}/\text{min}$
- (3) Temperature of a vaporizing room : 250 $^{\circ}\text{C}$
- (4) Temperature of a detector : 250 $^{\circ}\text{C}$
- (5) Carrier gas : N₂ 0.7 ml / min
- (6) Detector : FID
- (7) Range : 0.32 V
- (8) Sensitivity : 103 M Ω
- (9) Speed of chart : 5 mm / min

2.2 Gas chromatograph chart

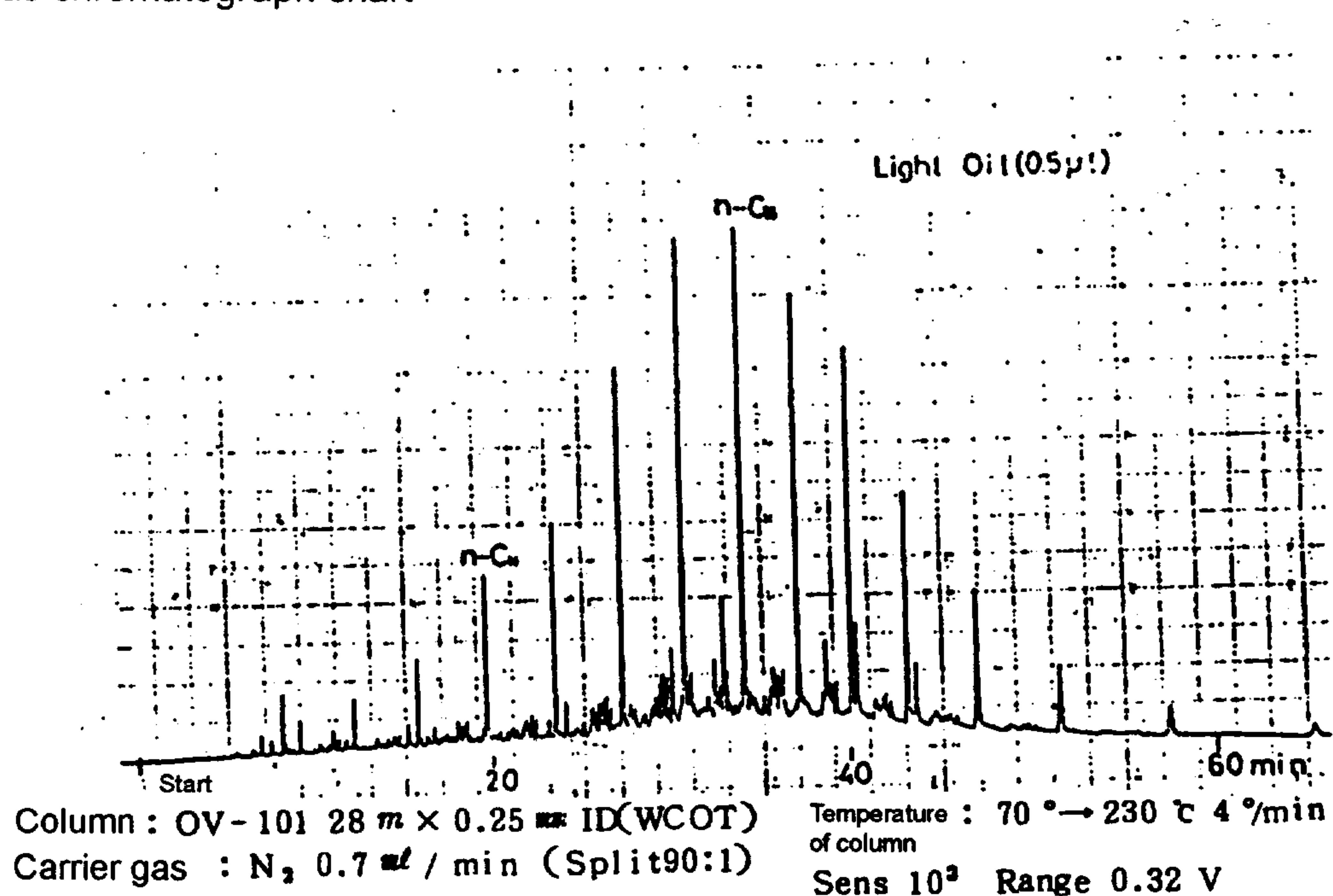


Figure 1. Gas chromatograph